
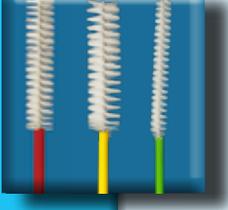
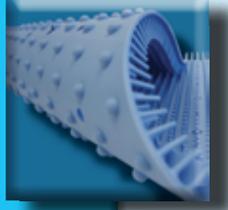


06/2020

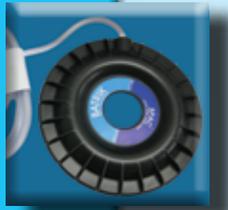

Catalogue 2020

What's New?

Valve Covers p. 23

FIRST CHECK
Beside Sponge kits
with Non-Enzymatic Detergent p. 45

Krystal View Ultrasound Gel p. 65

Interactive PDF Tip:

When on a Product page, click on these tabs to go back to the Table of Contents at any time!

Surgical Products & General Accessories

Suture Tags	04
Clamp Jaw Covers	05
Vascular Silicone Ties	06
Cautery Tip Cleaner	07
SKINS Markers	08
Pin Guards	09
Safe-Soak Mats & SPAC	10
Magnetic Mat	11
Tourniquets	12
Anti-Fatigue Mats	13
Stackable Setp Stools	17

Endoscopy Products

GOLFF - Anti-Fog Solution	20
Bite Blocks	21
Valve Covers	23
Channel Brushes	24
Valve/Control Brushes	27
Cleaning Brushes - V-Brush	28

Infection Control / Sterilization

Silicone Pin Mat	29
Corner Pockets	30
Instrument Protectors	31
STICS & ID Tapes	33
Tape Remover - Tape-B-Gone	34
Padlocks	35
Faceshields	39
Scrub Pads	40

NEW

Enzymatic Detergents	41
Instrument Lubricants	42
FIRST CHECK - Ready-To-Use Spray	43
FIRST CHECK - Sponge	44
FIRST CHECK - Beside Sponge kits with Non-Enzymatic Detergent	45
FIRST CHECK - Bedside Enzymatic Sponge Kits	46
FIRST CHECK - FIRST Wipes	47
FIRST CHECK - Transport Bags	48

Brushes

STERILIZATION BRUSHES

ARC Wire Brushes	49
Orthopedic Rasps Cleaning Brush	52
Stainless Steel Micro Brush	53
ARC Double Wire Brushes	54
Antimicrobial Twisted Wire Brushes	56

INSTRUMENT CLEANING BRUSHES

Flat Brushes	57
Precision Brush	58
Jar Brush	59
Sterilizer Drain Cleaning Brush	59
Acetabular Brush	60
All Purpose Brushes	61
Autoclave Cleaning Brush Kit	62


Women's Health

Hysterosonography Catheter	63
Ultrasound Gel	65




BATRIK Suture Clamp Tags are designed to protect and safely hold sutures without damage.

- Packed in a convenient foam base with easy peel-off (even with gloves) adhesive backing
- The foam base is affixed to a sterile area with the clamps presented in front of the Clamp Tags to allow for easy insertion
- Available **sterile** and **non sterile**
- Available in 2, 3 or 5 pairs
- Radio opaque
- Non-toxic
- Pre-counted
- Latex free

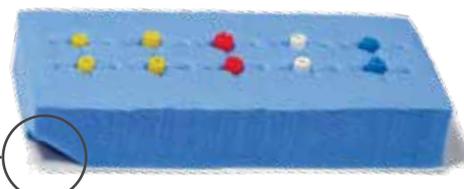

Easy Peel Corner

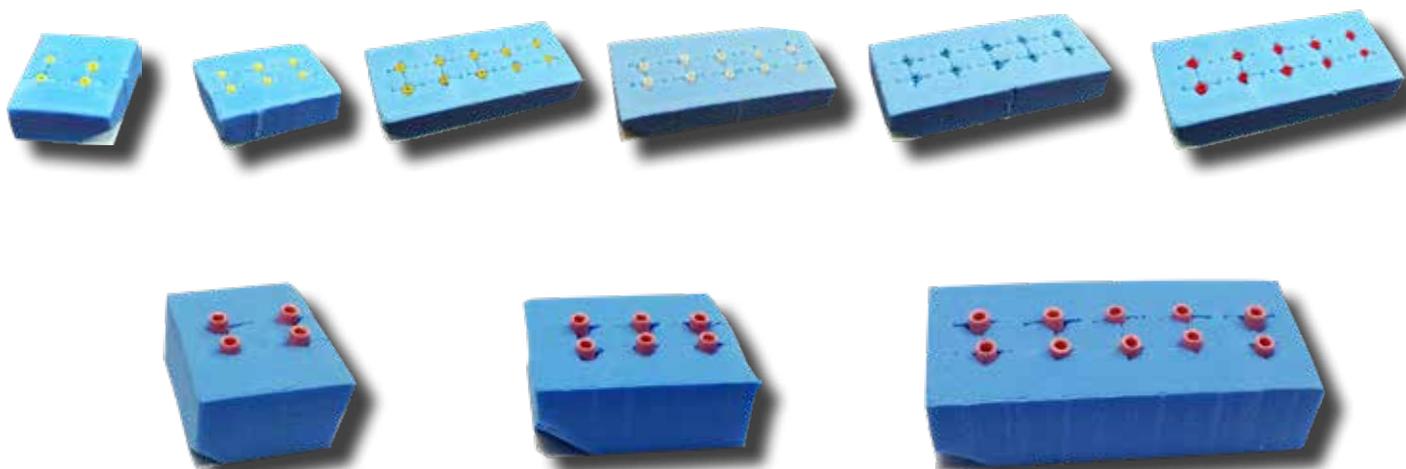

Suture Clamp Tags


Sterile item #	Non Sterile item #	Packaging	Size	Color
11-1002	911-1002	2 pairs / foam base	Standard ø 0.06" x 0.5" L	Yellow
11-1003	911-1003	3 pairs / foam base	Standard ø 0.06" x 0.5" L	Yellow
11-1005	911-1005	5 pairs / foam base	Standard ø 0.06" x 0.5" L	Yellow
11-1025	911-1025	5 pairs / foam base	Standard ø 0.06" x 0.5" L	Assorted
11-1045	911-1045	5 pairs / foam base	Standard ø 0.06" x 0.5" L	White
11-1065	911-1065	5 pairs / foam base	Standard ø 0.06" x 0.5" L	Blue
11-1085	911-1085	5 pairs / foam base	Standard ø 0.06" x 0.5" L	Red
11-3002	911-3002	2 pairs / foam base	Large ø 0.09" x 0.75" L	Pink
11-3003	911-3003	3 pairs / foam base	Large ø 0.09" x 0.75" L	Pink
11-3005	911-3005	5 pairs / foam base	Large ø 0.09" x 0.75" L	Pink

Our "Easy Peel Corner" is a feature unique to BATRIK

Packaging: Sterile: 1 foam base/pouch – 5 pouches/box
Non Sterile: sold per foam base

BATRIK Clamp Jaw Covers reduce trauma to vessels requiring occlusion, particularly when a hemostatic clamp is utilized. The thread-like indents allow for gentle, constant and non-slip tension.

- Manufactured from Medical-Grade silicone
- Gamma sterilized
- Available in 3 colors and 3 widths to fit most clamp sizes
- Convert almost any instrument to atraumatic
- Radio opaque
- Latex free

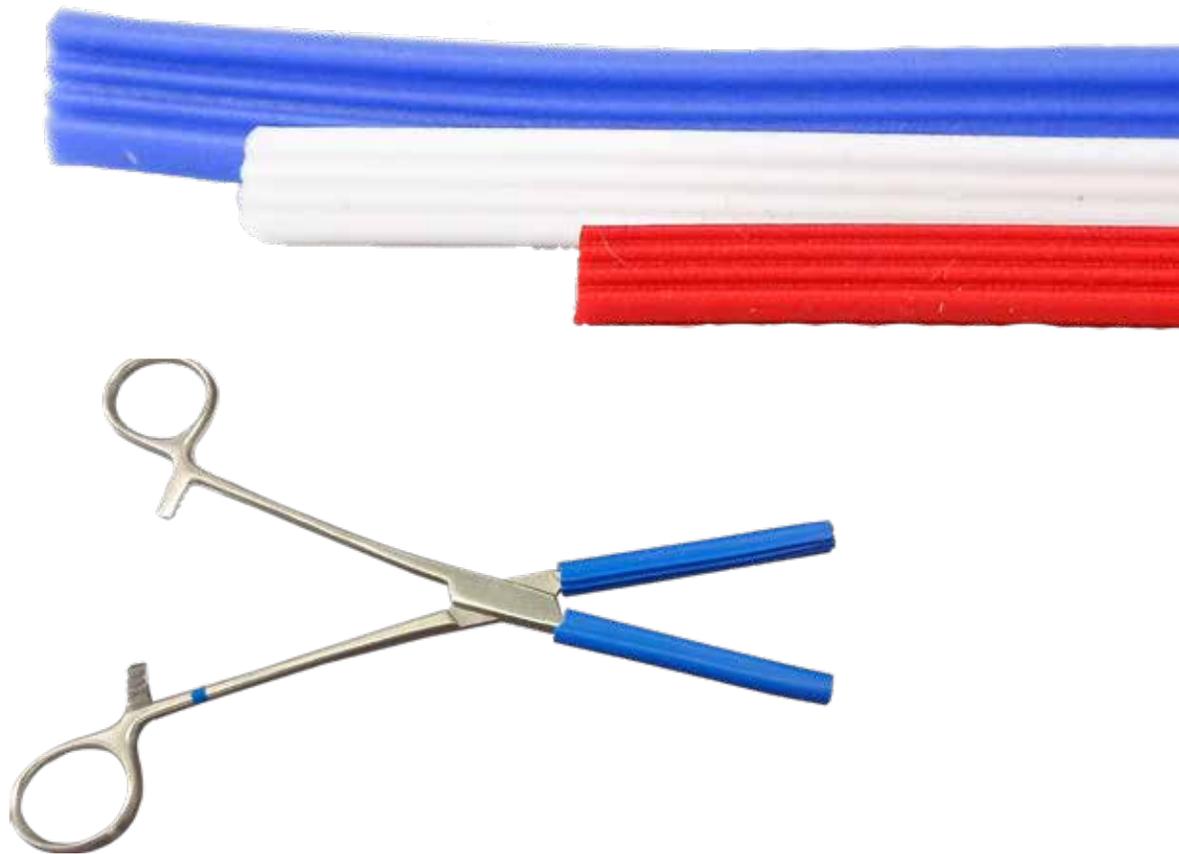

Clamp Jaw Covers

Sterile Item #	Non Sterile Item #	Dimensions – inches / cm	Color
12-1000	912-1000	0.1" x 5" / 0.25 x 12.7cm	Red
12-1001	912-1001	0.2" x 5" / 0.51 x 12.7cm	White
12-1002	912-1002	0.3" x 5" / 0.76 x 12.7cm	Blue

Packaging: Sterile: 2 clamps per sterile pouch – 10 pouches per box
Non Sterile: 20 clamps per bag (also available 100 units / bag)



BATRIK Vascular Silicone Ties permit safe and smooth retraction, occlusion and identification of vessels during surgical procedures.

Unlike clamps and retractors which can pinch the vessels flat, or even crack them, the BATRIK elliptical shape allows for the gentle handling of arteries, veins, nerves, tendons, spermatic cords, ureters, or whenever delicate retraction / occlusion is required.

- Manufactured of Medical-Grade silicone
- Do not rotate or twist on themselves when unpacked
- Easy to separate — eliminating the need for cardboard inserts
- Radio opaque
- Available in 4 colors and 4 sizes
- Latex free




Vascular Silicone Ties

Item #	Dimensions – inches / cm	Color
10-1000	0.075" x 16" / 0.2cm x 40cm	Mini-Red
10-1001	0.075" x 16" / 0.2cm x 40cm	Mini-Blue
10-1002	0.075" x 16" / 0.2cm x 40cm	Mini-Yellow
10-1003	0.075" x 16" / 0.2cm x 40cm	Mini-White
10-1010	0.120" x 16" / 0.3cm x 40cm	Maxi-Red
10-1011	0.120" x 16" / 0.3cm x 40cm	Maxi-Blue
10-1012	0.120" x 16" / 0.3cm x 40cm	Maxi-Yellow
10-1013	0.120" x 16" / 0.3cm x 40cm	Maxi-White
10-1101	0.040" x 16" / 0.1cm x 40cm	Extra Mini-Blue
10-1111	0.200" x 22" / 0.5cm x 56cm	Extra Maxi-Blue
10-1113	0.200" x 22" / 0.5cm x 56cm	Extra Maxi-White

Packaging: 2 ties per sterile pouch — 10 pouches per box

BATRIK Cautery Tip Cleaners assist in removing material from cautery tips for both mono-polar and bi-polar electrodes.

- Foam pad with integrated abrasive material
- Adheres to disposable and non-disposable drapes
- Only available non-sterile
- Radio opaque
- Latex free

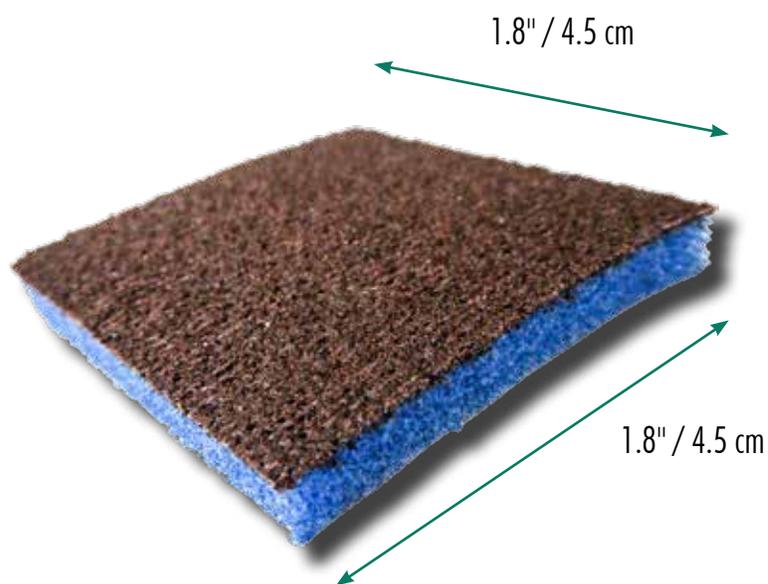

Cautery Tip Cleaners

Item #
60-10000

Description
Non Sterile - Cautery Tip Cleaners
1.8 x 1.8 in / 4.5 x 4.5 cm


Packaging
2200 / box

ver 04/04/2019



BATRIK SKINS™ – Surgical Skin Markers Designed with a beveled tip to enable the creation of multiple tracing lines (fine or bold) just by holding the pen at different angles.

- Made of lasting gentian violet ink that is non-toxic and non-tattooing
- Available in non-sterile for pre-operative use, or sterile for use in the operating room
- Available with or without 6 labels and a flexible 6" calibrated (in/mm) ruler
- Surgical prep resistant
- Latex free

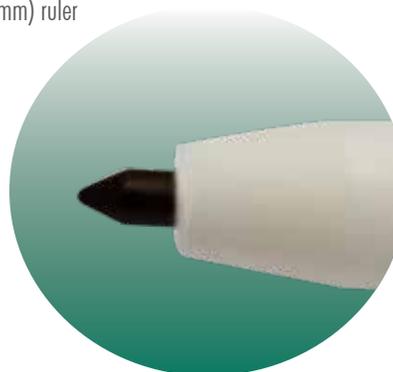

SKINS™ Surgical Skin Markers

Item #	Description	Packaging
S-3000	STERILE – Multi-Tip fine marker with ruler and label set	50 / box
S-3001	STERILE – Multi-Tip fine marker	50 / box
S-30005	STERILE – Multi-Tip fine marker with ruler and label set (Volume Packaging)	500 / box
S-30015	STERILE – Multi-Tip fine marker (Volume Packaging)	500 / box
S-93000	NON STERILE – Multi-Tip fine marker with ruler and label set	1000 / box
S-93001	NON STERILE – Multi-Tip fine marker	1000 / box
S-94000	NON-STERILE – Ruler	100 / bag
S-94001	NON-STERILE – Ruler (Bulk)	20,000 / box
S-95000	NON-STERILE – Set of 3 labels	100 / bag

SAVE by ordering our SKINS MARKERS VOLUME PACKAGING – more cost effective and reduces waste

ver 07/13/2016

BATRIK Pin Guards You will no longer have to use tape or tubing as our BATRIK Pin Guards will safely cover the sharp ends of your K-Wire and Steinmann Pins.

- Color coded by wire size which eliminates the selection of the wrong size
- A removable sleeve is provided to avoid injury and facilitate application (sterile only)
- Packaging designed to reduce waste
- Non sterile (100 units / bag) available with or without plastic sleeves
- Latex free



Pin Guards

Sterile item #	Non Sterile item # ⁽¹⁾	Color	Dimensions inches / mm	Sterile Packaging	Non-Sterile Packaging**
14-1028	914-1028	Yellow	0.028" / 0.71mm	80 units/box*	100 units/box
14-1035	914-1035	Light Blue	0.035" / 0.89mm	80 units/box*	100 units/box
14-1045	914-1045	White	0.045" / 1.14mm	80 units/box*	100 units/box
14-1054	914-1054	Dark Red	0.054" / 1.40mm	80 units/box*	100 units/box
14-1062	914-1062	Green	0.062" / 1.57mm	80 units/box*	100 units/box
14-1078	914-1078	Pink	0.078" / 1.98mm	80 units/box*	100 units/box
14-1093	914-1093	Grey	0.093" / 2.36mm	80 units/box*	100 units/box
14-1109	914-1109	Beige	0.109" / 2.77mm	80 units/box*	100 units/box
14-1125	914-1125	Purple	0.125" / 3.17mm	80 units/box*	100 units/box
14-1028-1	N/A	Yellow	0.028" / 0.71mm	40 units/box	N/A
14-1035-1	N/A	Light Blue	0.035" / 0.89mm	40 units/box	N/A
14-1045-1	N/A	White	0.045" / 1.14mm	40 units/box	N/A
14-1062-1	N/A	Green	0.062" / 1.57mm	40 units/box	N/A
14-1141	914-1141	Brown	0.141" / 3.58mm	80 units/box*	100 units/box
14-1156	914-1156	Dark Blue	0.156" / 3.96mm	80 units/box*	100 units/box
14-1188	914-1188	Red	0.188" / 4.77mm	80 units/box*	100 units/box

*Sterile Packaging: 2 units per sterile pouch; 40 pouches per box.

** Non-Sterile: Does not include plastic sleeve

⁽¹⁾ Non-sterile sold without plastic sleeves. See price list for part number and price with sleeve.

ver 12/30/2016



BATRIK Safe-Soak™ Absorbent Floor Mats are designed to absorb fluids from the operating room floor quickly and efficiently, to ensure that the floor beneath remains slip-free, safe and dry during surgical procedures. A skid-resistant backing is available with the Blue Absorbent Mat. All mats are single use and non-sterile.

- Absorbs a variety of fluids
- Provides quick and easy clean-up
- Highly absorbent
- Latex free, disposable, and non-sterile
- Conforms to uneven surfaces



Safe-Soak™ Absorbent Floor Mats

Item #	Description	Size	Packaging	Absorbency
SS-20-1000	Safe-Soak™ / White extra absorbent mat	28" x 48" (71cm x 122cm)	20 per box	4.5 liters / mat
SS-20-2000	Safe-Soak™ / Blue absorbent mat with skid-resistant backer	28" x 40" (71cm x 102cm)	20 per box	2.5 liters / mat
SS-20-3000	Safe-Soak™ / Yellow absorbent mat	30" x 52" (76cm x 132cm)	15 per box	3.0 liters / mat

BATRIK Safe-Soak™ SPAC floor suction disk, quickly removes fluid waste from the operating room floor through the application of suction. These suction disks have a highly absorbent base material which allows quick fluid wicking and removal when suction is applied. Safe-Soak™ suction disks are latex free and non sterile.

- Easily maneuverable around operating room floor
- fluid is captured from bottom through a highly absorbent material
- opening permits fluid collection that come into contact with its surface
- highly absorbent and fast wicking
- latex free, disposable, and non-sterile
- highly durable and won't be damaged if stepped on
- the opening allows for easy movement using your foot
- the hose is securely affixed on its side to reduce risk of catching equipment



Safe-Soak™ SPAC


Item #	Description	Size	Packaging
SS-21-1000	Safe-Soak™ SPAC	8.625" / 22cm diameter	12 / box

ver 06/22/2018

BATRIK Magnetic Mats provide a **sterile, neutral field** for the safe placement and **hands-free transfer** of instruments.

- **Raised magnets** hold instruments firmly in place while allowing easy grasping of instruments
- **Round** magnets to avoid silicone rubber perforations
- **Flexible** and **non-slip** backing drape conforms to any surface shape

Reusable and Autoclavable (up to 134°C)
Latex free

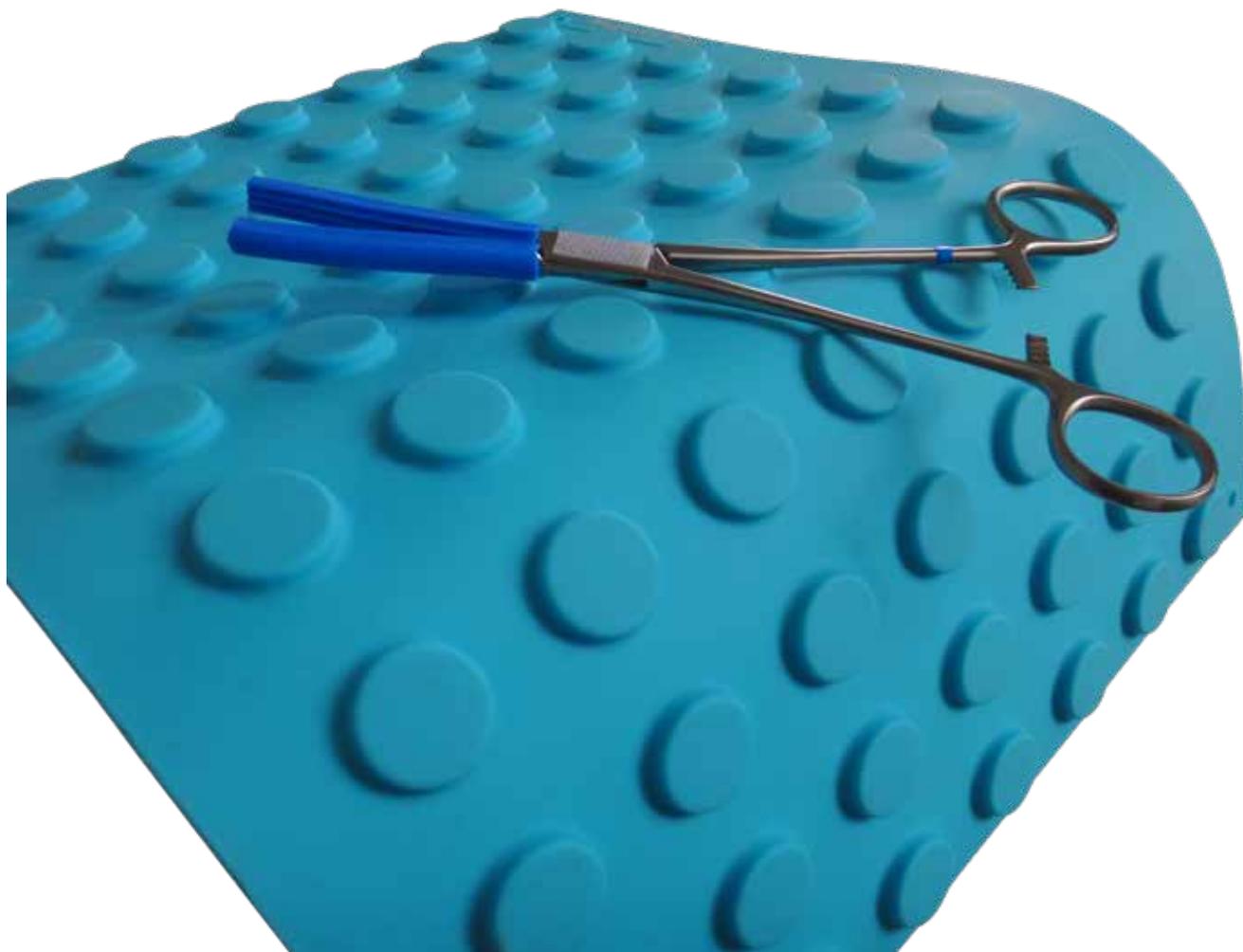

ver 07/13/2016

Magnetic Mats

Item #	Dimensions – inches / cm	Packaging
MAG-02-1	12" x 15" / 30cm x 38cm	1 unit
MAG-02-5	12" x 15" / 30cm x 38cm	5 units



BATRIK Tourniquets are comfortable and easy to use on all patients, including; the elderly, undernourished, children and newborns.

- Safeguards hospital staff and patients from latex allergens
- **New design:** Our material will feel like latex, but has a more comfortable texture for the patients and less stress on the skin
- Available in different sizes
- Compact sizes is convenient for most phlebotomy trays
- Latex free

Single Use Offers the Following Benefits :

- most tourniquets are impossible to disinfect and re-usable units are often contaminated
- will reduce the risk of cross contamination (patient to patient or patient to healthcare provider)
- may help reduce the spread of nosocomial infections
- may improve your facility’s compliance with OSHA’s 1910-1030 (d)(1) “Blood-borne Pathogens Standard”
- 5-6% of the North American population have a latex allergy, particularly children with spina bifida



Facilitate single-use policy with our **25 units per roll**, which comes in a convenient dispenser box for easy access!

TOURNY™ + Tourniquets

Item #	Description	Color	Dimensions	Packaging
K-002	TOURNY™ + Non Sterile	Blue	1" x 18" / 2.5cm x 45cm	1 unit / bag
K-004	TOURNY™ + Non Sterile	Blue	1" x 18" / 2.5cm x 45cm	100 units / box
K-009	TOURNY™ + Non Sterile	Blue	0.200" x 22" / 0.5cm x 56cm	20 units / box
K-010	TOURNY™ + Non Sterile	White	0.200" x 22" / 0.5cm x 56cm	20 units / box
K-400	TOURNY™ + Non Sterile	Blue	1" x 18" / 2.5cm x 45cm	16 rolls - 400 units/box
K-400PB	TOURNY™ + Non Sterile (Pebbled)	Light Blue	1" x 18" / 2.5cm x 45cm	16 rolls - 400 units/box
K-100	TOURNY™ + Non Sterile	Assorted	1" x 18" / 2.5cm x 45cm (3 rolls of each Blue, Orange, Green and Pink)	12 rolls - 300 units/box
K-254	TOURNY™ + Non Sterile (Pebbled)	Light Blue	1" x 18" / 2.5cm x 45cm	250 units/box
K-2-18RL	TOURNY™ + Non Sterile	Blue	1" x 18" / 2.5cm x 45cm	1 roll
K-2-18PBRL	TOURNY™ + Non Sterile (Pebbled)	Light Blue	1" x 18" / 2.5cm x 45cm	1 roll
K-3-18RL	TOURNY™ + Non Sterile	Green	1" x 18" / 2.5cm x 45cm	1 roll
K-6-18RL	TOURNY™ + Non Sterile	Pink	1" x 18" / 2.5cm x 45cm	1 roll
K-7-18RL	TOURNY™ + Non Sterile	Orange	1" x 18" / 2.5cm x 45cm	1 roll
K-2-29RL	Long TOURNY™ + Non Sterile	Blue	1" x 29.5" / 2.5cm x 75cm	1 roll


ver 06/20/2019



BATRIK Medical Blue Cloud-Gel™ Anti-Fatigue Mat is truly the medical industry standard. It has been scientifically engineered using patented gel technology combined with premium foam to provide the ultimate in comfort and support to surgeons, nurses and technicians.

Both the fabrication and material composition is specifically designed for use in all hospital applications, particularly where movement is restricted for prolonged durations. Our material composition will withstand the repeated exposure to a wide range of fluids.

Feel the **Blue Cloud-Gel™** Advantage designed for the O.R. !

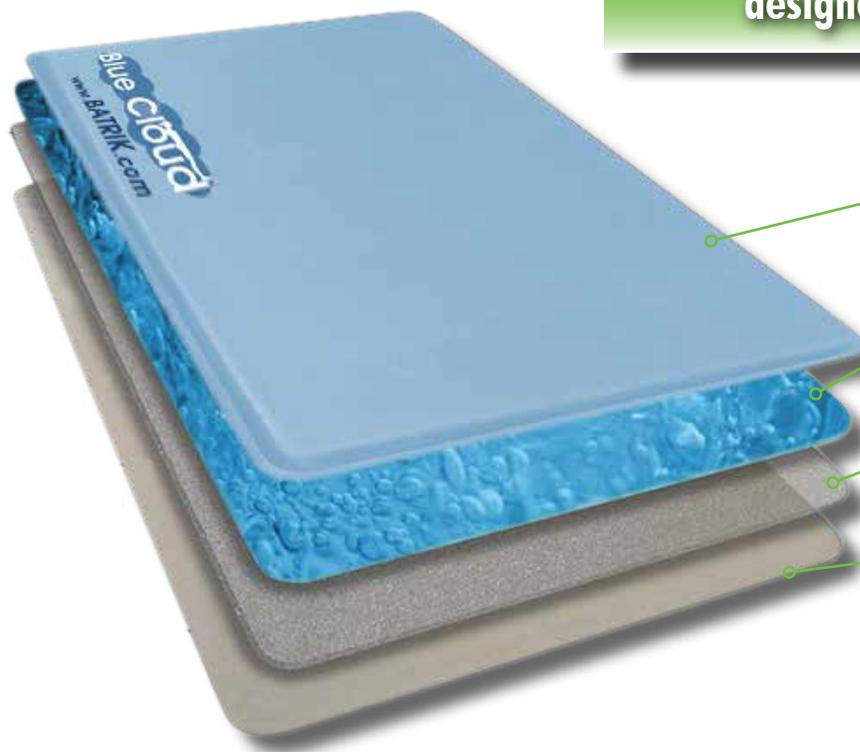

Both top and bottom are stain resistant / anti-microbial & anti-bacterial

Integrated **Blue Cloud-Gel™**

Superior Quality Foam

Anti-Slip Backing

 Latex-Free

Cart Washer Compatible (Contact Us for more details)

COMPARE THE BENEFITS AND FEEL THE DIFFERENCE.

- Integrated **Blue Cloud-Gel™** will not “Bottom-Out” like conventional foam mats.
- Properly conforms and supports feet while promoting efficient circulation, reducing physical discomfort and fatigue.
- Promotes proper posture while diminishing back, knee and foot discomfort.
- Improves comfort and ergonomics while also enhancing energy.
- Fabricated with the high standard for the operating room, but the anti-fatigue properties will benefit every department.
- Studies have demonstrated that anti-fatigue mats improve performance, energy and post-operative discomfort.

www.BlueCloudGel.com

ver 06/20/2019





FEATURES INCLUDE:

- Unique integrated **Blue Cloud-Gel™** technology provides superior long lasting anti-fatigue advantages over conventional foam based products.
- Incorporated high performance foam.
- Both the top and bottom are stain resistant and anti-microbial / anti-bacterial.
- The bottom material is Non-Slip Certified by the National Floor Safety Institute.
- Proprietary 3/4" Batrik **Blue Cloud-Gel™** Foam technology.
- Excellent friction properties when damp or dry surgical booties are worn.
- Water-resistant and non-absorbent.
- Easy to clean with most hospital quaternary cleaners.*
- US standards: CAL 117 and NFPA 260 Compliant.
- Reusable
- 1 Year Warranty



Blue Cloud-Gel Anti-Fatigue Mat

Item #	Dimensions – inches / cm	Packaging
800-2401	18 x 24 / 46 x 61	1 / Box
800-2405	18 x 24 / 46 x 61	5 / Box
800-3201	20 x 32 / 51 x 81	1 / Box
800-3203	20 x 32 / 51 x 81	3 / Box
800-4801	20 x 48 / 51 x 122	1 / Box
800-4802	20 x 48 / 51 x 122	2 / Box
800-7201	30 x 72 / 76 x 183 (special order)	1 / Box


*With the wide array of cleaners, we recommend that you test your solution beforehand.

Did You Know? We also offer the excellent anti-fatigue mat "Nimbus", using only our superior foam based material (NOT for O.R. Application). You will be amazed with the benefits and our comfort technology.

www.BlueCloudGel.com



BATRIK Medical's Blue Cloud-Nimbus™ Anti-Fatigue Mat proprietary high-density "energy-return" polyurethane materials will provide unmatched comfort, relief and ergonomic support.

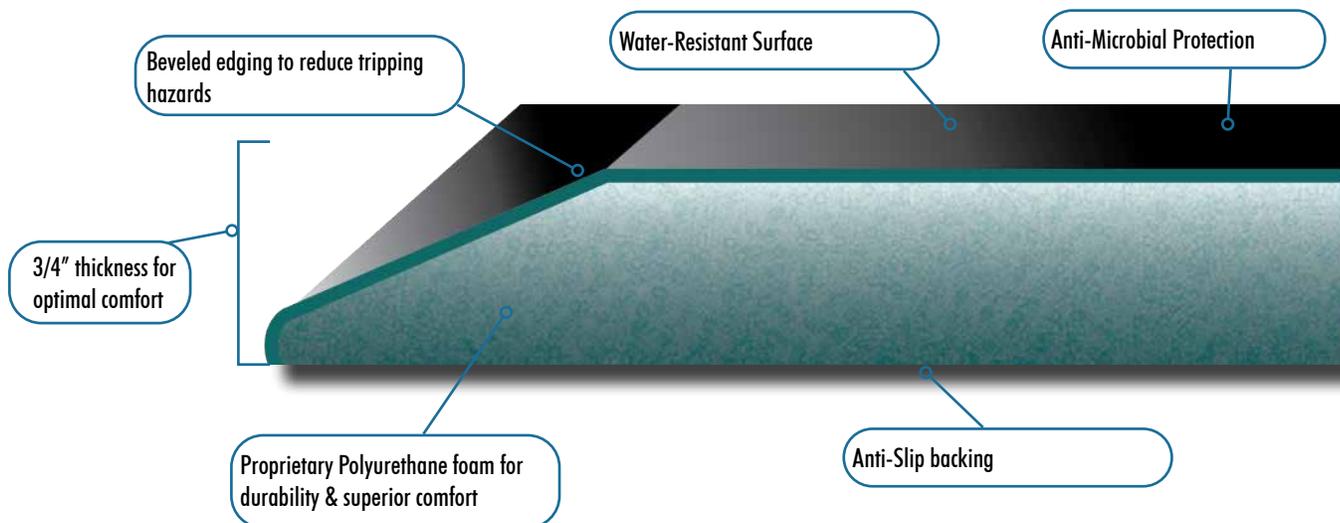

COMPARE THE BENEFITS AND SEE THE DIFFERENCE

- Our Energy-Return properties will help reduce leg, foot, lower back discomfort and work-related injuries.
- Feel the difference compared to any other foam or rubber based anti-fatigue mat on the market as our bounce-back memory restores energy.
- Longer lasting support as 3/4" (1.9 cm) proprietary foam blend will not bottom out like most other foam based anti-fatigue mats.
- Eco-Friendly materials are used, so the **Blue Cloud-Nimbus™** is not just good for you, but also good for the environment.
- If within 30 days you are not 100% satisfied with the added benefits, contact your dealer for a refund.*



ver 06/20/2019

www.BlueCloudGel.com



FEATURES INCLUDE

- Anti-Microbial properties blended into the mat to provide added protection.
- Highly Water-Resistant Surface will not absorb water or dirt. (Not commonly found on other typical foam based mats.)
- One piece construction ensures longer useful life.
- Textured design on surface and high-traction base ensures greater safety. (US Standard NFSI 101-C)
- Durable design is ideal for commercial applications, yet still ergonomically & aesthetically pleasing for home or office environments.
- Flammability/Burn resistant per CFR 16 Parts 1630/1631
- Non-Toxic and Latex Free
- Light weight and easy to clean.**
- 5 Year Warranty.***



Blue Cloud Nimbus Anti-Fatigue Mat

Item #	Dimensions – inches / cm	Packaging	Color	Item #	Dimensions – inches / cm	Packaging	Color
809-3201	20 x 32 / 51 x 81	1 / Box	Black	810-3201	20 x 32 / 51 x 81	1 / Box	Taupe
809-3204	20 x 32 / 51 x 81	4 / Box	Black	810-3204	20 x 32 / 51 x 81	4 / Box	Taupe
809-4801	20 x 48 / 51 x 122	1 / Box	Black	810-4801	20 x 48 / 51 x 122	1 / Box	Taupe
809-4804	20 x 48 / 51 x 122	4 / Box	Black	810-4804	20 x 48 / 51 x 122	4 / Box	Taupe
809-3601	24 x 36 / 61 x 91	1 / Box	Black	810-3601	24 x 36 / 61 x 91	1 / Box	Taupe
809-7201	20 x 72 / 51 x 183	1 / Box	Black	810-7201	20 x 72 / 51 x 183	1 / Box	Taupe
809-6001	36 x 60 / 91 x 152	1 / Box	Black	810-6001	36 x 60 / 91 x 152	1 / Box	Taupe

NEW PRODUCTS

- 809-1705 Reusable Anti-Fatigue Mat for our step stool OR-36363 (p.24)¹
- 899-1710 Single-Use Anti-Fatigue Mat for our step stool OR-36363 (p.24)¹

¹With adhesive strip to secure it in place.

* Refund applies only to the purchase price, excludes any transport, handling or taxes.
 ** We recommend using most soap based products. Do not submerge the mat in water.
 *** 5 Year Warranty for Non-Industrial usage; 4 Year Warranty for Industrial usage.

ver 06/20/2019



Stackable Step Stools by IMPERIAL SURGICAL

BATRIK is proud to be distributing on behalf of IMPERIAL SURGICAL the unique, antimicrobial Stackable Step Stool. Made of Aluminium, it has truly been designed with Infection Prevention as a prime focus. It is Cart Washable and can even be sterilized.*



The IMPERIAL SURGICAL Stackable Step Stools are designed for the O.R., but are equally useful hospital wide. Made of light weight aluminum, our Stackable Step Stools are protected with **an anti-microbial powder coating**. Furthermore, they can support more weight than most of our competitors (see Figure 1).

- **Improved** model
- **Anti-slip surface** for better stability
- **No rubber corners:** this was a source of previous fluid collection (see the picture on next page, circled in red)
- **No welds** or large hollow opening to trap debris or fluids
- **Transport cart** Available (Figure 2)
- **NEW:** You may now purchase a fully autoclavable stepstool without the mat! See last page for order details.



Color-code
Your Stepstools
by
Department!

* The OR-36364 can be sterilized. See last page for details.

Can support the weight of half a Jeep Wrangler!



Figure 1



Figure 2

Transport Cart
(see last page for order details)

ver 05/05/2020

www.StackingStepStool.com



www.surgmed.com • info@surgmed.com • +1 514-631-7988 • Fax +1 514-631-9083

CAN THE COMPETITION COMPARE?

Anti-Microbial Coating


Patent pending technology. These stools are often in contact with moisture, fluid and blood; and will benefit from anti-microbial protection. The benefits include:

- Reduced risks of infections¹
- Prevents the spread of dangerous micro-organisms and suppresses the growth of mold, bacteria, yeast and mildew
- Provides superior durability and protection for essentially the life of the product.²

No Corner Rubbers

Did you know that the corner rubbers found on some of the competition's stackable stepstools can become a source for bacteria? With time, they erode and create deep cavity channels that can collect fluids and thus, create an environment for bacteria to propagate.

IMPERIAL SURGICAL has removed these rubber corners for your safety.

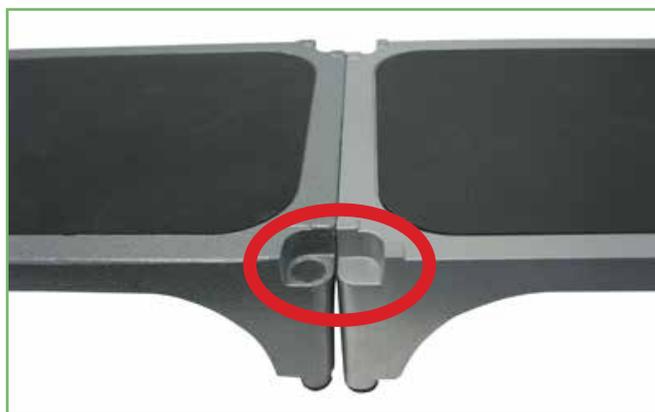

Compatible

If you have older models or similar looking competitor stepstools, these newly designed models will stack and inter-lock perfectly.

High density foot pegs

Our foot pegs are made of denser material than the competition, and will not break or come loose easily. And if they do; we will replace them free of charge! Only Transport fees will apply (see Warranty).

Now Available

Reusable and Single-use (O.R.)* Anti-Fatigue Mats. They provide comfort, support and reduce leg & back fatigue when standing for extended period.

*Single-Use Replacement Mats to reduce cross contamination when needed.

7 years warranty

We offer one of the best warranty in the market. Our 7 years warranty covers all parts (step stool only).

World-Class Customer Service

Your satisfaction is our utmost priority.

ver 06/20/2019

ORDERING INFORMATION

Material: Cast Aluminium

Weight: 8.3 lbs (3.8 kg)

Weight Capacity: 300 lbs (136 kg)

Dimensions: 17.4" x 13" x 4.7" (L x W x H)

Packaging (1): 19" x 14" x 6" (10 lbs)

Packaging (4): 26.5" x 20" x 15" (45 lbs)



OR-36363

Stacking Step Stools

Item #	Description
OR-36363	Stackable Stepstool with mat
OR-36363-02	Stackable Stepstool with mat
OR-36363-03	Stackable Stepstool with mat
OR-36363-07	Stackable Stepstool with mat
OR-36364	Stackable Stepstool without mat

Color

- Black
- Dark Blue
- Green
- Orange

Autoclavable



OR-3693

Optional Accessory

Item #	Description	Packaging
OR-3693	Transport Cart	1
809-1705	Reusable Anti-Fatigue Mat*	5 units / box
899-1710	Single-Use Anti-Fatigue Mat*	10 units / box

*With adhesive strips to secure it in place.



899-1710

Replacement Parts

Item #	Description
OR-36363-1A	Replacement Pad for OR-36363
OR-36363-1B	Replacement Foot Pegs

ver 06/20/2019

Visit us at www.StackingStepStool.com or
 Write us at info@surgmed.com
 for more product details

BATRIK GOLFF™ Anti Fog Solution Laparoscopic surgical procedures often have problems with fogging. The BATRIK Sterile Anti-fog Solution, GOLFF™, is the solution that sets the standard for all other products.

- Ensures clearer vision and avoids unnecessary interruptions during surgery
- Prevents the occurrence of fog on lenses
- Adhesive radiopaque sponge for secure placement
- Easy squeeze 6 mL bottle
- Latex free
- OEM (available)



GOLFF™


Item #	Description	Packaging
GOLFF	Adhesive sponge and anti-fog per set	20 sterile sets / box

ver 06/26/2019

BATRIK Endoscopy BO2™ Bite Blocks with Nasal Oxygen Delivery

significantly enhances patient safety by providing supplemental oxygenation during gastrointestinal endoscopy preventing hypoxia in sedated patients, while ensuring superior patient comfort.

- Lower procedural costs by obviating the need for separate nasal prongs and oxygen masks, the BO2 Bite Block helps save time and money
- Ensure patient comfort and safety as with other Bite Blocks from BATRIK Medical, the BO2 Bite Block is designed to be free of sharp edges to ensure patient comfort and safety. Additionally, our bite blocks are guaranteed to be latex-free, lead-free and DEHP-free.
- Choose from disposable and reusable options. The BO2 Bite Block is available in single use and reusable options. Use our flexible ARC Wire Brush #41-32505, which will allow you to properly clean the nasal prongs and inlet of the reusable BO2.
- Ensure secure fit. An optional retention strap is available to ensure a secure fit. These straps, like the bite block, are latex free, lead free and DEHP free.

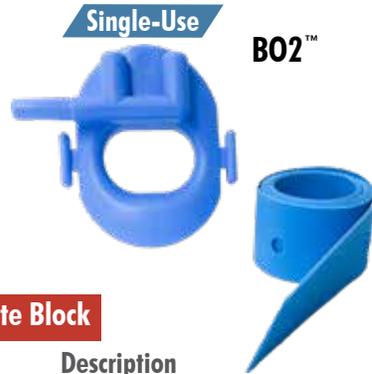

BO2™ Bite Block

Item #	Description	Packaging
BO2	18mm / 54FR - Single use, includes strap	100 / Box
BO2-N	18mm / 54FR - Single use, without strap	100 / Box
BO2-RN	18mm / 54FR - Reusable, without strap	25 / Box
BO2-R	18mm / 54FR - Reusable, includes strap	25 / Box

BATRIK Endoscopy Bitestrong™ Bite Blocks

are designed to provide a safe passage through the patient's mouth during TEE thus protecting the patient's teeth as well as the endoscope. Strict quality control procedures and the careful selection of raw materials used in the manufacture of our Bite Blocks is your guarantee that they are latex free, lead free and DEHP free.



Bite Blocks

Item #	Description	Packaging
B4	19mm / 60FR - Single use	100 / Box
B4T	19mm / 60FR - Single use, includes strap	100 / Box
B4T-N	19mm / 60FR - Single use, without strap	100 / Box
B1-25	16mm / 48FR - Pediatric single use, includes strap	25 / Box
B1	16mm / 48FR - Pediatric single use, includes strap	100 / Box
B1-N	16mm / 48FR - Pediatric single use, without strap	100 / Box



BATRIK Endoscopy Bitestrong™ Single Use Bite Blocks

BEST SELLER

- Our molding process ensures a smooth texture finish with no sharp edges (flashing) that can often irritate a patient’s mouth
- Our bite blocks are individually wrapped in sealed, easy-to-open PE bags
- Available without dental retention
- Additional straps are also available, if required




Bitestrong™ Single Use Adult Bite Block

Item #	Description	Packaging
B3	20mm / 60FR - Single use, includes strap	100 / Box
B3-N	20mm / 60FR Reusable, without strap	100 / Box
B3-SA	Same as B3, with strap <u>attached</u> to one side	100 / Box
B3-NDR	Same as B3, no dental retention rim	100 / Box
B3-NDR-N	Same as B3-N, no dental retention rim, without strap	100 / Box

BATRIK Endoscopy Bitestrong™ Reusable Bite Blocks

- Our reusable Bite Blocks are designed for repeated usage and will withstand commonly practiced chemical and steam sterilization cycles (134°C / 275°F)⁽¹⁾
- Our molding process ensures a smooth textured finish with no sharp edges (flashing) that can often irritate a patient’s mouth
- Our Bite Blocks are individually wrapped in sealed, easy-to-open PE bags
- Available with or without an adjustable retention strap
- Additional straps are also available, if required

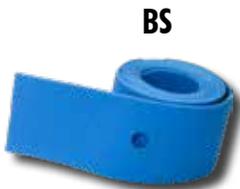

Bitestrong™ Reusable Adult Bite Block and Straps

Item #	Description	Packaging
B6	19mm / 60FR Reusable	100 / Box
B6T	19mm / 60FR Reusable, includes strap	100 / Box
B6T-N	19mm / 60FR Reusable, without strap	100 / Box
B6-25	19mm / 60FR Reusable	25 / Box
B6T-25	19mm / 60FR Reusable, includes strap	25 / Box
B6T-N-25	19mm / 60FR Reusable, without strap	25 / Box
BS	Single-Use Bite-Block straps 55cm x 2cm / 21.7" x 0.79"	100 / Box

Order your single use Bite Block sample box today! **Item # B-SAMPLE-S**
The BATRIK BITESTRONG BiteBlock Single Use Sample Box includes:

Order your reusable Bite Block sample box today! **Item # B-SAMPLE-R**
The BATRIK BITESTRONG BiteBlock Reusable Sample Box includes:

ver05/05/2020



BATRIK Couvrette™ Valve Covers eliminate the need for manual cleaning and reprocessing of reusable valves. Our Single-Use Valve Covers effectively reduce cross contamination between patients and also reduce fluid leakage from the parts during procedures.

- Easy and secure fit to all valves
- Reduce the potential for errors
- Single-use
- Non sterile
- Latex free



Reusable valve covers are difficult to clean. Save time and reduce risk.

Biopsy Valve Covers

Item #	Description	Color	Packaging
16-0022	Biopsy Valve Cover - Olympus® & Fujinon® (G5 Series et Les plus recentes)	Blue	100 / Box
16-0044	Biopsy Valve Cover - Pentax®	Red	100 / Box

ver 05/05/2020



BATRIK ENDOSCOPY Reusable Channel Brushes DAYBRUSH™


If a facility still prefers a Reusable brush, then the DayBrush™ is the perfect solution.

BATRIK has developed a disposable channel brush that can be used for cleaning endoscopes for an entire day. You no longer need to worry if your present reusable brush is cleaning effectively. After processing ten (10) endoscopes⁽¹⁾, simply discard the brush.

The unique structure of the Reusable DayBrush™ enables the following

- Soft tip ensures easy manipulation within the channel with no risk of damage
- The back-end is sealed and softened to prevent fluid from getting trapped inside the tube. The softening process also prevents damage if the cleaning brush is used in reverse
- Can be steam sterilized
- The volume of bristles on the DayBrush™ is 20mm in length — almost double most other brushes. A large brush provides more surface contact during the cleaning process
- The BATRIK Reusable DayBrush™ has been tested to effectively clean up to 10 scopes⁽¹⁾
- Significant per use savings when compared to reusable or disposable channel cleaning brushes

⁽¹⁾ Discard after cleaning 10 scopes or sooner if any brush degradation is noticed. After ten applications the brush MAY have the appearance that it can still effectively clean, but we discourage prolonged use due to risk of cross contamination. Actual number of uses will vary based on handling and processing.

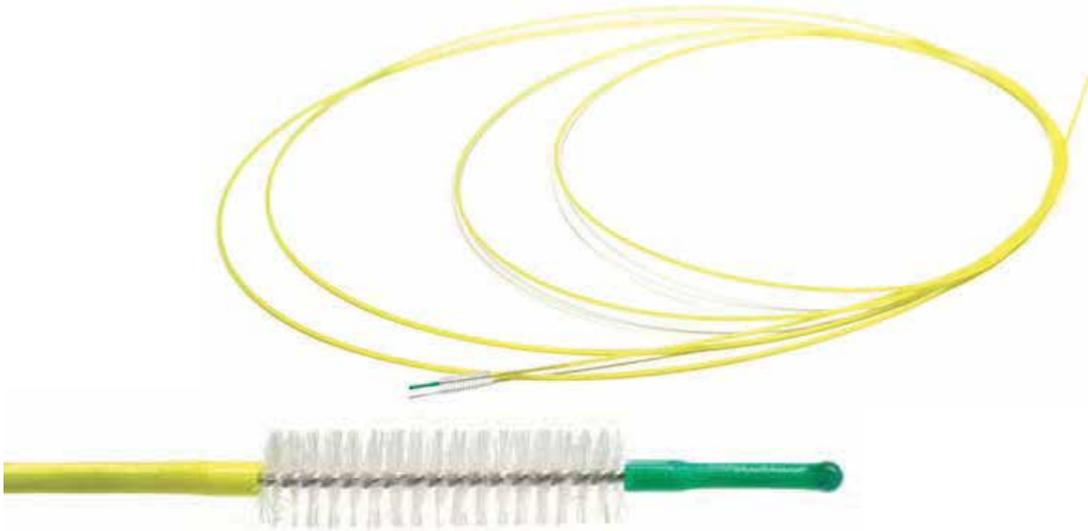

ADD THE SUFFIX...

F = receive a kit with our Flat Sponge (Packaging: Box of 40). Applicable on certain models only, call for details.

K = receive a kit with the single use Valve/Control Cleaning Brush 5/11 Double Brushes (item # 27-57787). Applicable on all models.

FK = includes both sponge and valve cleaning brush.



5mm Cylindrical DayBrush™

Item #	Length inches / cm	Brush Length inches / mm	Brush (Ø) inches / mm	Sheath (Ø) inches / mm	Brush Ends
30-18050-T	70" / 180cm	0.79" / 20mm	0.20" / 5mm	0.07" / 1.75mm	Single
30-23050-T	90" / 230cm	0.79" / 20mm	0.20" / 5mm	0.07" / 1.75mm	Single
30-23050-2T	90" / 230cm	0.79" / 20mm	0.20" / 5mm	0.07" / 1.75mm	Double

Packaging: 50 / Box

ver 02/01/2017


BATRIK ENDOSCOPY Single Use ARC™ Channel Cleaning Brushes

BATRIK offers one of the largest selection of Endoscopy Channel Cleaning Brushes in the industry. Our brushes are available in a variety of lengths, material, and styles to accommodate the widest range of endoscopic devices used in the healthcare industry today.

No matter what your requirements, you'll be able to find the brush that suits your specific needs. If you don't see the brush you require, just call us. We have well over 50 different channel brushes available and if you can't find the exact one that you need, no problem — we can custom fabricate to your specifications.

Why Choose a BATRIK ARC Channel Brush?

- One of the industry's longest brush heads on the market: 15mm to 24 mm in length. More bristles translate into more cleaning contact within the channels
- Bristle fan tip on the head of the brush protects and cleans the channel. No sharp end or hard metal tips that can scratch or cause damage
- Latex free
- Available in two different materials in order to accommodate different customer preferences:
 - POLYPROPYLENE is the industry standard as it is a low cost material and offers greater flexibility.
 - DELRIN is slightly more expensive, but due to its more rigid feeling, it is easier to insert and manipulate. This in turn reduces time and labor which can offset the added cost. ALL BATRIK small diameter brushes are supplied only with DELRIN, otherwise they would be too fragile

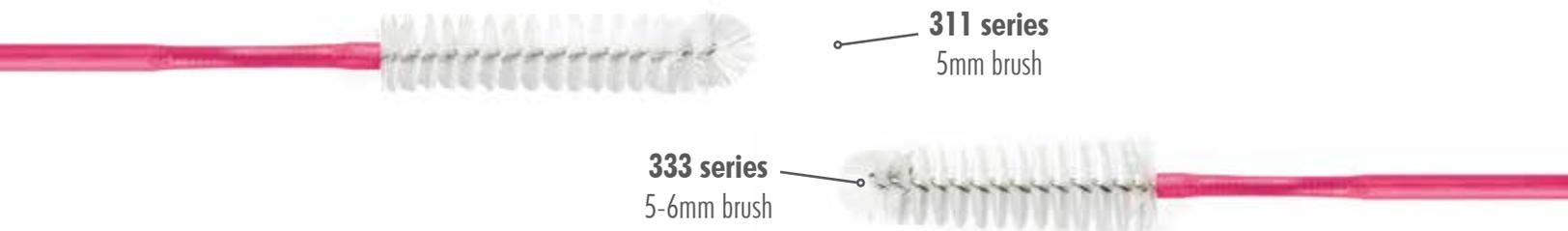

F= receive a kit with our Flat Sponge (Packaging: Box of 40). Applicable on certain models only, call for details.

K= receive a kit with the single use Valve / Control Cleaning Brush 5/11 Double Brushes (item # 27-57787). Applicable on all models.

FK= includes both sponge and valve cleaning brush.

5 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.85 MM – PUNCHY PINK POLYPROPYLENE

Item #	Length inches / cm	Brush Length inches / mm	Brush (Ø) inches / mm	Sheath (Ø) inches / mm	Brush Ends
311-23050	90" / 230cm	0.94" / 24mm	0.20" / 5mm	0.073" / 1.85mm	Single
311-23050-2	90" / 230cm	0.94" / 24mm	0.20" / 5mm	0.073" / 1.85mm	Double
311-18050	70" / 180cm	0.94" / 24mm	0.20" / 5mm	0.073" / 1.85mm	Single
311-10050	39" / 100cm	0.94" / 24mm	0.20" / 5mm	0.073" / 1.85mm	Single

Packaging: 50 / Box

Please note that you may receive the brush in a RED color instead of PINK. It is perfectly normal as we are slowly phasing out the RED color. The quality of the product remains identical (31 Series).

5/6 MM TAPERED BRUSHES – SHEATH Ø 1.85 MM – PUNCHY PINK POLYPROPYLENE


Item #	Length inches / cm	Brush Length inches / mm	Brush (Ø) inches / mm	Sheath (Ø) inches / mm	Brush Ends
333-23056	90" / 230cm	0.79" / 20mm	0.20, 0.24" / 5, 6mm	0.073" / 1.85mm	Single
333-23056-2	90" / 230cm	0.79" / 20mm	0.20, 0.24" / 5, 6mm	0.073" / 1.85mm	Double

Packaging: 50 / Box

Please note that you may receive the brush in a RED color instead of PINK. It is perfectly normal as we are slowly phasing out the RED color. The quality of the product remains identical (33 Series).

ver 08/23/2018

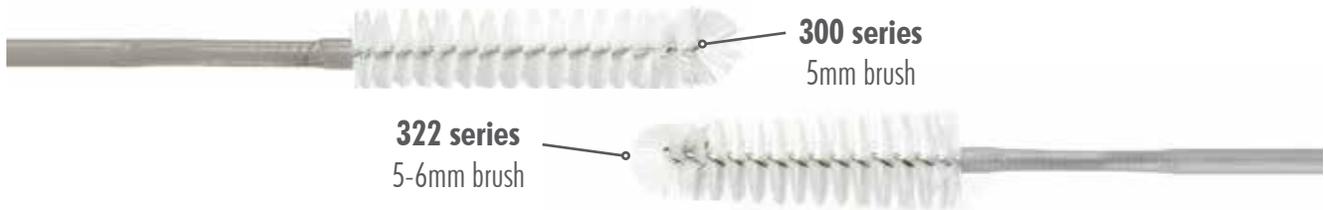

BATRIK ENDOSCOPY Single Use ARC™ Channel Cleaning Brushes**5 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.75 MM – METALLIC GREY DELRIN**

Item #	Length inches / cm	Brush Length inches / mm	Brush (Ø) inches / mm	Sheath (Ø) inches / mm	Brush Ends
300-18050	70" / 180cm	0.79" / 20mm	0.20" / 5mm	0.070" / 1.75mm	Single
300-23050	90" / 230cm	0.79" / 20mm	0.20" / 5mm	0.070" / 1.75mm	Single
300-23050-2	90" / 230cm	0.79" / 20mm	0.20" / 5mm	0.070" / 1.75mm	Double

Packaging: 50 / BoxPlease note that you may receive the brush in a **WHITE** color instead of **GREY**. It is perfectly normal as we are slowly phasing out the **WHITE** color. The quality of the product remains identical (30 Series)**5/6 MM TAPERED BRUSHES – SHEATH Ø 1.75 MM – METALLIC GREY DELRIN**

Item #	Length inches / cm	Brush Length inches / mm	Brush (Ø) inches / mm	Sheath (Ø) inches / mm	Brush Ends
322-23056	90" / 230cm	0.79" / 20mm	0.20, 0.24" / 5, 6mm	0.070" / 1.75mm	Single
322-23056-2	90" / 230cm	0.79" / 20mm	0.20, 0.24" / 5, 6mm	0.070" / 1.75mm	Double

Packaging: 50 / BoxPlease note that you may receive the brush in a **WHITE** color instead of **GREY**. It is perfectly normal as we are slowly phasing out the **WHITE** color. The quality of the product remains identical (32 Series)**3 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.55 MM – METALLIC PURPLE DELRIN**

Item #	Length inches / cm	Brush Length inches / mm	Brush (Ø) inches / mm	Sheath (Ø) inches / mm	Brush Ends
344-23030	90" / 230cm	0.79" / 20mm	0.12" / 3mm	0.060" / 1.55mm	Single
344-18030	70" / 180cm	0.79" / 20mm	0.12" / 3mm	0.060" / 1.55mm	Single
344-10030	39" / 100cm	0.79" / 20mm	0.12" / 3mm	0.060" / 1.55mm	Single

Packaging: 50 / BoxPlease note that you may receive the brush in a **GREEN** color instead of **PURPLE**. It is perfectly normal as we are slowly phasing out the **GREEN** color. The quality of the product remains identical (34 Series)**2.5 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.00 MM – YELLOW DELRIN**

Item #	Length inches / cm	Brush Length inches / mm	Brush (Ø) inches / mm	Sheath (Ø) inches / mm	Brush Ends
35-23025	90" / 230cm	0.59" / 15mm	0.098" / 2.5mm	0.039" / 1.00mm	Single

Packaging: 50 / Box**2 MM CYLINDRICAL BRUSHES – SHEATH Ø 1.00 MM – YELLOW DELRIN**

Item #	Length inches / cm	Brush Length inches / mm	Brush (Ø) inches / mm	Sheath (Ø) inches / mm	Brush Ends
35-23020	90" / 230cm	0.59" / 15mm	0.079" / 2.0mm	0.039" / 1.00mm	Single

Packaging: 50 / Box

ver 08/23/2018

BATRIK Valve / Control Cleaning Brushes

VALVE / CONTROL CLEANING BRUSH – 5/7



Item #	Length inches / cm	Brush Length inches / mm	Brush 1 (Ø) inches / mm	Brush 2 (Ø) inches / mm	Description
27-05197	5.00" / 13 cm	1.00" / 25 mm	0.20" / 5 mm	0.28" / 7 mm	Dble. Trumpet

Packaging: (Single-use recommended) Bag of 3 - **ITEM #** 27-05197-50 / Bag of 50

VALVE / CONTROL CLEANING BRUSH – 5/11



Item #	Length inches / cm	Brush 1 Length inches / mm	Brush 2 Length inches / mm	Brush 1 (Ø) inches / mm	Brush 2 (Ø)x inches / mm	Description
27-57787	6.00" / 15cm	1.18" / 30 mm	0.79" / 20 mm	0.43" / 11 mm	0.20" / 5 mm	Dble. Universal

Packaging: (Single-use recommended) Bag of 3 - **ITEM #** 27-57787-50 / Bag of 50

VALVE / CONTROL CLEANING REUSABLE BRUSH – 5/11 – EXTRA LONG



Item #	Length inches / cm	Brush 1 Length inches / mm	Brush 2 Length inches / mm	Brush 1 (Ø) inches / mm	Brush 2 (Ø) inches / mm	Description
40-18511	8.07" / 20.5cm	1.57" / 40 mm	0.79" / 20 mm	0.43" / 11 mm	0.20" / 5 mm	Dble. Universal

Packaging: (Re-usable and can be steam sterilized) Bag of 3

VALVE / CONTROL CLEANING BRUSH – SINGLE USE BRUSH



or



Item #	Length inches / cm	Brush Length inches / mm	Brush (Ø) inches / mm	Description
36-10010*	4.3" / 11 cm	1.57" / 40 mm	0.39" / 10 mm	Single end large

Packaging: (Single-use recommended) Bag of 3 - **ITEM #** 36-10010-50 / Bag of 50

*Actual brush handle may vary, but brush head is unchanged.

ver 05/22/2018

BATRIK ENDOSCOPY V-Brush - Reusable External Cleaning Brush External cleaning brush for endoscopes / optical instruments. With its unique shape and innovative design, the BATRIK Reusable External Cleaning Brush is the ultimate tool for cleaning the exterior of your endoscopic instruments without the risk of damaging expensive equipment.

- Easy-grip handle ensures comfort and control
- Numerous dense bristles provide additional surface contact
- Superior endoscopic scope cleaning with a brush specifically designed for this purpose



Handle Length / Diameter
120mm / 12mm

Brush Length / Diameter
90mm / 25mm

Total Length
210mm

- Polyamide (TYNEX®) bristles
- Superior bristle durability
- Steamable thermoplastic handle
- Ergonomic handle
- Steam / chemical sterilizable
- Latex free

Reusable External Cleaning Brush

Item #	Description	Packaging
36-20010	Reusable External Cleaning Brush	3 units / Bag

BATRIK ENDOSCOPY V-Brush Single Use External Cleaning Brush External cleaning brush for endoscopes / optical instruments. An ideal fit for most endoscopic instruments, the BATRIK Single Use External Cleaning Brush will provide excellent cleaning capabilities. It's U-shape bristle design provides continuous and effective cleaning. There are no sharp or metallic ends that can scratch or cause damage to valuable equipment.

Many of the same features as the reusable brush

- Polypropylene bristles
- Vinyl caps
- Ergonomic handle
- Latex free

Handle Length / Diameter
80mm / 3mm



Brush Length / Diameter
70mm / 25mm

Total Length
150mm

Single Use External Cleaning Brush

Item #	Description	Packaging
36-20005	Single Use External Cleaning Brush	50 units / Box



ver 06/29/2018

BATRIK Silicone Pin Mat cushions, protects and organizes delicate surgical instruments for sterilization. Residual moisture collection below the silicone mat is significantly reduced due to our unique **elevated "Nub"** design. These "Nubs" allow for greater air-flow ventilation and drying.

New Features

- Now with **more holes** to facilitate better venting and reduce water collection
- Unique **elevated nubs** also provide superior venting and a reduction in water collection
- **Custom sizes** available on request (call for details)

Reusable and Autoclavable
Latex free

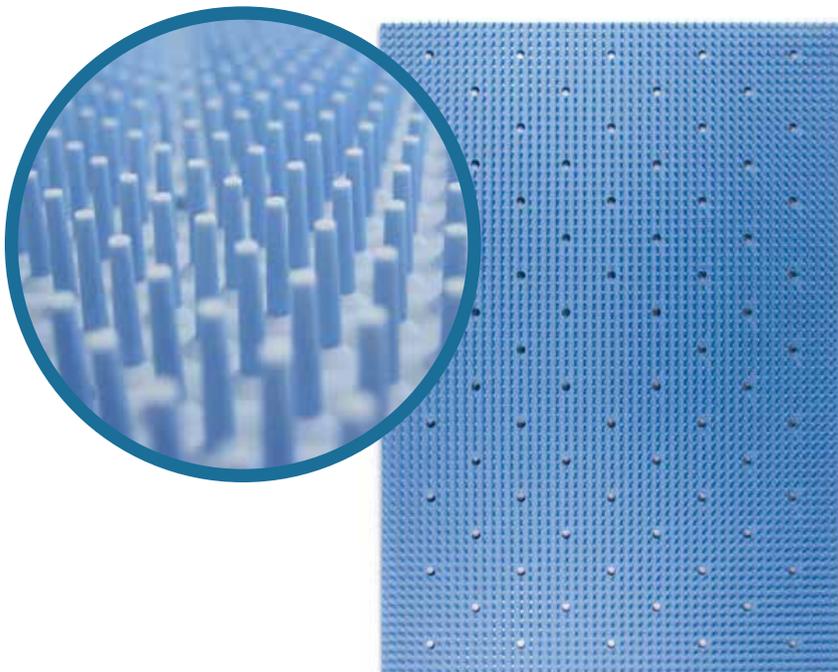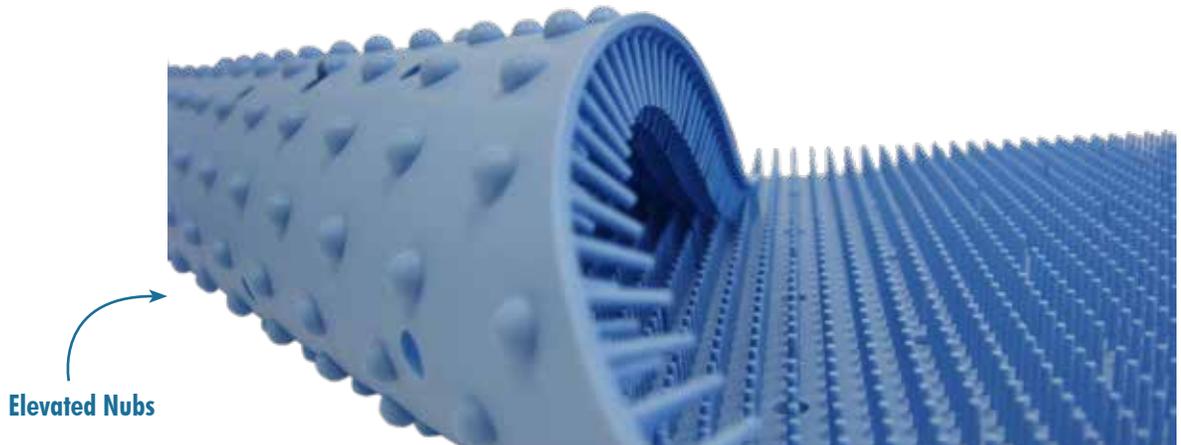

ver 06/29/2018

Silicone Pin Mat

Item #	Dimensions – inches / cm	Packaging
13-2001	13" x 18" / 33cm x 45cm	1 unit
13-2005	13" x 18" / 33cm x 45cm	5 units



BATRIK Corner Pockets™ Avoid frustration and costs from having to reprocess your surgical trays due to ripped sterile wrapping. Corner Pockets™ will prevent tearing associated with sharp corners or feet on surgical trays. They will fit all types of corners (even the round ones), and the bright orange-red color makes them easily identifiable and will reduce the risk of being discarded. We are also pleased to introduce one of the only non-foam single-use corner protectors (see below).

- No moisture retention
- Reusable or single-use options
- Non sterile
- Latex free



The Reusable corner pockets can withstand repeated sterilization up to 134°C; after 10 cycles, no product degradation was observed.

Corner Pockets™

Item #	Description	Packaging
18-CP07	Corner Pockets™ — Reusable	80 pcs / Bag
18-CORNER	Corner Pockets™ — Reusable	4 Bags / Box (320 pcs)
18-EVAL7	Corner Pockets™ — Reusable	24 pcs / Bag
18-CP02*	Corner Pockets™ — Single-Use	240 pcs / Bag
18-EVAL2*	Corner Pockets™ — Single-Use	24 pcs / Bag



AN INDUSTRY FIRST – SINGLE USE CORNER POCKETS!

Our Single-Use corner pockets have the following advantages:

- save time; **easy** to use and position
- **no moisture** retention

** Although the product appears similar to the reusable device, please be advised that it is blended differently to reduce costs and has only been validated as single use. Re-using the single use Corner Pocket may compromise sterility.*

***NEW* FOAM CORNER PROTECTORS**

Soft Pockets™

Item #	Description	Packaging
18-CPF-02	Soft Pockets™ - Single-Use Foam Corner Protectors	80 sets / Box

“Our foam is 0.25” thicker than the competitors. Especially useful with heavier Instrument sets. Our size difference is hardly noticeable, but the benefits are very apparent.”
(24 boxes per pallet)



6” x 6” x 1.50”
(15.2cm x 15.2cm x 3.8cm)

ver 06/28/2018




BATRIK Instrument Protectors will protect instrument tips and edges from damage during handling and sterilization.

- Available in a wide array of sizes for almost every application
- Sizes are color coded
- Opaque and tinted protectors are radiopaque
- Steam and gas permeable
- Non sterile
- Latex free



Variations Include

- Perforated protectors (holes) which reduce vacuum seal and moisture build-up
- Osteotome protectors
- Double protectors
- Retractor protectors
 - "Weitlaner" Perforated "Square" Style
 - Single-Tip Perforated Clear "Tube" Style
- Validated for pre-vacuum sterilization

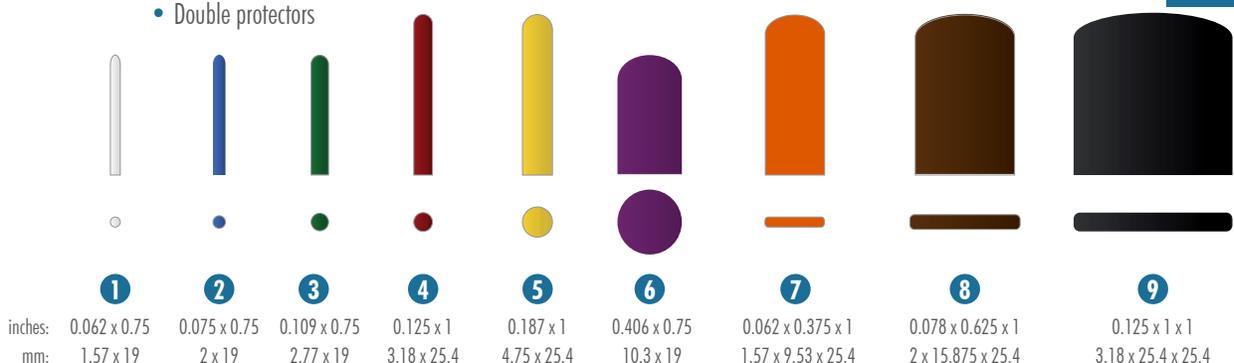

Opaque – Solid Color

Item # (regular)	Item # (perforated)
1 RS-01	PS-01
2 RS-02	PS-02
3 RS-03	PS-03
4 RS-04	PS-04
5 RS-05	PS-05
6 RS-06	PS-06
7 RS-07	PS-07
8 RS-08	PS-08
9 RS-09	PS-09
RS-10	PS-10

Packaging: 100 units/Bag

Tinted Color

Item # (regular)	Item # (perforated)	Color
RT-01	PT-01	White
RT-02	PT-02	Blue
RT-03	PT-03	Green
RT-04	PT-04	Red
RT-05	PT-05	Yellow
RT-06	PT-06	Maroon
RT-07	PT-07	Orange
RT-08	PT-08	Brown
RT-09	PT-09	Black
RT-10	PT-10	Assorted

Packaging: 100 units/Bag

Clear - No Color

Item # (regular)	Item # (perforated)
1 RC-01	PC-01
2 RC-02	PC-02
3 RC-03	PC-03
4 RC-04	PC-04
5 RC-05	PC-05
6 RC-06	PC-06
7 RC-07	PC-07
8 RC-08	PC-08
9 RC-09	PC-09
RC-10	PC-10

Packaging: 100 units/Bag



ver 07/13/2016




Double Protectors

Item #	Dimensions – inches / mm	Color
DG-01	0.0625 x 0.75" / 1.59 x 19.05mm	White
DG-02	0.078 x 0.75" / 1.98 x 19.05mm	Blue
DG-03	0.109 x 0.75" / 2.78 x 19.05mm	Green
DG-04	0.125 x 1" / 3.18 x 25.4mm	Red
DG-05	0.1875 x 1" / 4.76 x 25.4mm	Yellow
DG-07	0.0625 x 0.375 x 1" / 1.59 x 9.53 x 25.4mm	Orange
DG-10	Assorted	Assorted



Packaging: 50 units/Bag - Perforated only

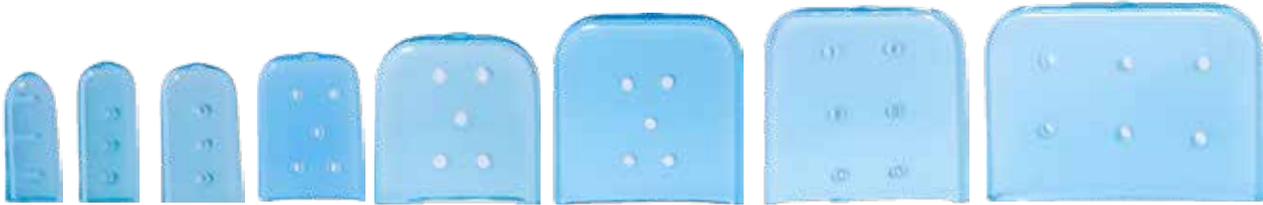


Osteotome Protectors

Item #	Dimensions – inches / mm
OP-11	0.25 x 0.75" / 6.35 x 19.05mm
OP-12	0.375 x 1" / 9.53 x 25.4mm
OP-13	0.5 x 1" / 12.7 x 25.4mm
OP-14	0.75 x 1" / 19.05 x 25.4mm
OP-15	1 x 1" / 25.4 x 25.4mm
OP-16	1.25 x 1.25" / 31.75 x 31.75mm
OP-17	1.5 x 1.5" / 38.1 x 38.1mm
OP-18	1.75 x 1.25" / 44.45 x 31.75mm
OP-20	Assorted

Retractor Protectors Single-Tip

Item #	Dimensions – inches / cm	Packaging
RGS-01	1.5" / 3.8cm	55 units / Bag
RGS-02	2" / 5.1cm	45 units / Bag
RGS-03	3" / 7.62cm	40 units / Bag

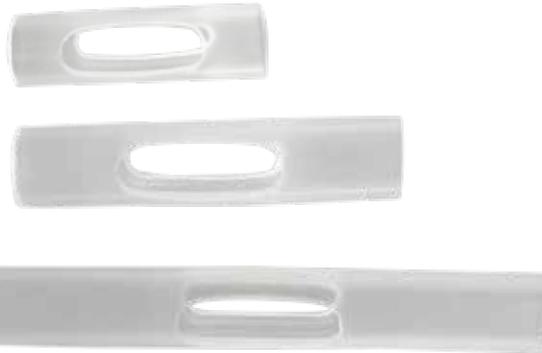

Packaging: 50 units/Bag - Perforated only

Retractor Protectors "Weitlaner"

Item #	Size	Color
RGM-01	Small	Pink
RGM-02	Small	Clear
RGM-03	Medium	Green
RGM-04	Medium	Clear



Packaging: 50 units/Bag

ver 07/13/2016

BATRIK S.T.I.C.S.™ a Simple Tracking Identification Color System

- Quickly identifies instruments by department or procedure
- Easy to apply
- Long lasting
- Tape adheres securely
- A free protective cover is provide with every sheet
- Latex free
- Non sterile

Each 8.5 x 11" (216 x 279mm) sheet contains the following adhesive strips:

- 1" (25.4mm) x 220
- 1.5" (38.1mm) x 44
- 2" (50.8mm) x 44
- total tape per sheet 374" (9500mm)

Solid Colors

Item #	Color	
T-850-2-01	White	
T-850-2-02	Royal Blue	
T-850-2-03	Kelly Green	
T-850-2-04	Red	
T-850-2-05	Yellow	
T-850-2-06	Purple	
T-850-2-07	Orange	
T-850-2-08	Brown	
T-850-2-09	Black	
T-850-2-10	Gray	
T-850-2-11	Lavender	
T-850-2-12	Beige	
T-850-2-13	Light Blue	
T-850-2-14	Lime Green	
T-850-2-15	Pink	
T-850-2-17	Plum	
T-850-2-18	Mint Green	

Black Stripes

Item #	Color	
T-850-3-01	Black / White	
T-850-3-03	Black / Kelly Green	
T-850-3-04	Black / Red	
T-850-3-05	Black / Yellow	
T-850-3-07	Black / Orange	
T-850-3-10	Black / Gray	
T-850-3-11	Black / Lavender	
T-850-3-12	Black / Beige	
T-850-3-13	Black / Light Blue	
T-850-3-14	Black / Lime Green	
T-850-3-15	Black / Pink	
T-850-3-17	Black / Plum	

White Stripes

Item #	Color	
T-850-4-02	White / Royal Blue	
T-850-4-03	White / Kelly Green	
T-850-4-04	White / Red	
T-850-4-05	White / Yellow	
T-850-4-07	White / Orange	
T-850-4-08	White / Brown	
T-850-4-11	White / Lavender	
T-850-4-13	White / Light Blue	
T-850-4-14	White / Lime Green	
T-850-4-15	White / Pink	
T-850-4-17	White / Plum	

Tri-Color Stripe

Item #	Color	
T-850-5-51	Red / White / Blue	
T-850-5-54	Red / Black / Green	

Specialty Stripes

Item #	Color	
T-850-7-04	Yellow / Red	

Four Color

Item #	Color
T-850-2-31	Red / Yellow / Purple / Orange
T-850-2-32	Blue / White / Black / Green

BATRIK Identification Tapes are a simple and inexpensive identification solution for marking surgical instruments.

- Quickly identifies instruments by department or procedure
- Easy to apply
- Withstands repeated exposure to sterilization
- Has more tape per roll (300" / 762cm)
- Are available in 9 colors and 2 widths (0.25" / 0.635cm and 0.125" / 0.3175cm)
- Packaged in a convenient clear box
- Non sterile

*sterilization reports are available. All samples tested validated the efficiency of sterilization, even underneath the tape

Identification Tapes

Item #	Item #	Color
0.125" / 0.3175 cm	0.25" / 0.635 cm	
T-125-01	T-250-01	White
T-125-02	T-250-02	Blue
T-125-03	T-250-03	Green
T-125-04	T-250-04	Red
T-125-05	T-250-05	Yellow
T-125-06	T-250-06	Purple
T-125-07	T-250-07	Orange
T-125-08	T-250-08	Brown
T-125-09	T-250-09	Black



ver 07/13/2016



BATRIK Color Coding System Padlocks Security Tags

Plastic seals for instrument containers and general hospital applications are specifically designed for the healthcare market. Our material selection, fabrication process and engineering ensures consistent quality and high temperature compliance with the various steam sterilization temperatures in world wide hospital applications.




- **NEW:** Now available in different sizes!
- Large base for easier handling
- Easytwist breaking design
- Fully open design (45°) facilitates insertion into instrument containers
- Will not melt or soften with any established steam sterilization protocol
- Made of polypropylene
- Latex and DEHP Free
- Custom Colors and OEM available upon request




Simply fold the base and pull to break the padlock

Now available in 3 different formats


ACTUAL SIZE

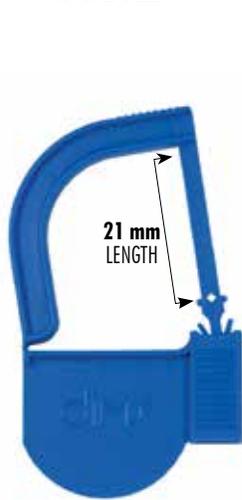

700 Series
"Original"

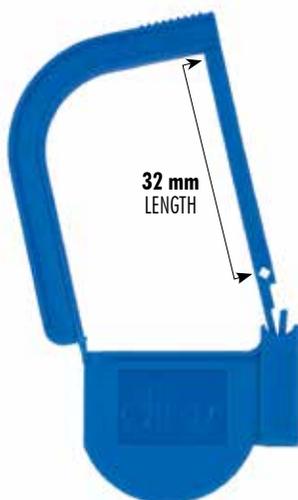

800 Series
"NEW"

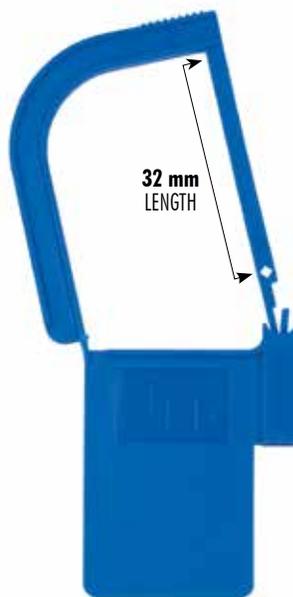

900 Series
"NEW"



With Indicator Dot



Printed Padlocks



Numbered

800 & 900 Series have a single looking hook for use with smaller diameters

ver 08/29/2018



BATRIK Color Coding System Padlocks Security Tags

The 700 series padlocks is our original design, perfect for “Aesculap” type of containers. They also work for containers with a short, large opening.



Double locking hooks can be inserted in a large guide hole

The grooves on the padlock facilitates handling



700-Series Padlocks

Regular padlocks



Color	100 Units	200 Units	500 units
White	15-701	15-711	15-710
Blue	15-702	15-712	15-720
Green	15-703	15-713	15-730
Red	15-704	15-714	15-740
Yellow	15-705	15-715	15-750
Pink	15-706	15-716	15-760
Orange	15-707	15-717	15-770
Black	15-709	15-719	15-790
Light Blue	15-722	15-732	15-742

With Indicator dots



Color	100 Units	200 Units	500 units
White	15-701-IND	15-711-IND	15-710-IND
Blue	15-702-IND	15-712-IND	15-720-IND
Green	15-703-IND	15-713-IND	15-730-IND
Red	15-704-IND	15-714-IND	15-740-IND
Yellow	15-705-IND	15-715-IND	15-750-IND
Pink	15-706-IND	15-716-IND	15-760-IND
Orange	15-707-IND	15-717-IND	15-770-IND
Black	15-709-IND	15-719-IND	15-790-IND
Light Blue	15-722-IND	15-732-IND	15-742-IND

Printed padlocks




Color	100 Units
White “CLEAN”	15-701C
Blue “STERILE”	15-702ST
Red “REPAIR”	15-704R
Yellow “SHARPEN”	15-705S
Yellow “LOANER”	15-705L
Black “DISCARD”	15-709D

Numbered Padlocks



We also offer Custom printed padlocks upon request (up to 8 characters).
Call us for more details!

ver 03/07/2019

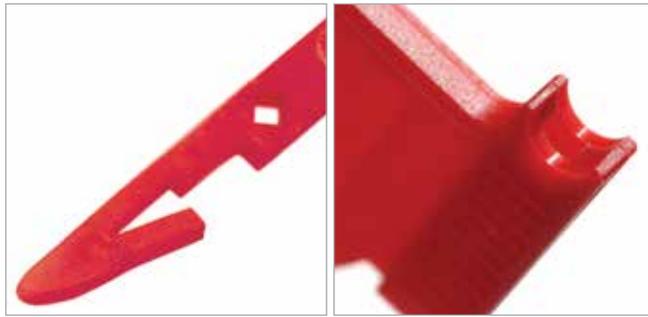

Simple and fine locking hook can be inserted in a large guide hole

The locking system allows for easy closure

BATRIK Color Coding System Padlocks Security Tags

The 800 series padlocks have a longer arm (32 mm) and have a single, thinner hook than the 700 series. They are perfect for almost any type of containers: even those with long and thin channels.



800-Series Padlocks

Regular padlocks



Color	100 Units	200 Units	500 units
White	15-801	15-811	15-810
Blue	15-802	15-812	15-820
Green	15-803	15-813	15-830
Red	15-804	15-814	15-840
Yellow	15-805	15-815	15-850
Pink	15-806	15-816	15-860
Orange	15-807	15-817	15-870
Black	15-809	15-819	15-890
Light Blue	15-822	15-832	15-842

With Indicator dots



Color	100 Units	200 Units	500 units
White	15-801-IND	15-811-IND	15-810-IND
Blue	15-802-IND	15-812-IND	15-820-IND
Green	15-803-IND	15-813-IND	15-830-IND
Red	15-804-IND	15-814-IND	15-840-IND
Yellow	15-805-IND	15-815-IND	15-850-IND
Pink	15-806-IND	15-816-IND	15-860-IND
Orange	15-807-IND	15-817-IND	15-870-IND
Black	15-809-IND	15-819-IND	15-890-IND
Light Blue	15-822-IND	15-832-IND	15-842-IND

ver 03/07/2019

BATRIK Color Coding System Padlocks Security Tags

The 900 series padlocks are identical to the 800 series, but with a larger rectangular base, which allows for the placement of a 20x30mm tracking labels.

- Ideal to add labels/Indicators or markings to the extra large base



The locking system allows for easy closure



Large Rectangular base



900-Series Padlocks

Regular padlocks



Color	100 Units	200 Units	500 units
White	15-901	15-911	15-910
Blue	15-902	15-912	15-920
Green	15-903	15-913	15-930
Red	15-904	15-914	15-940
Yellow	15-905	15-915	15-950
Pink	15-906	15-916	15-960
Orange	15-907	15-917	15-970
Black	15-909	15-919	15-990
Light Blue	15-922	15-932	15-942

With Indicator dots



Color	100 Units	200 Units	500 units
White	15-901-IND	15-911-IND	15-910-IND
Blue	15-902-IND	15-912-IND	15-920-IND
Green	15-903-IND	15-913-IND	15-930-IND
Red	15-904-IND	15-914-IND	15-940-IND
Yellow	15-905-IND	15-915-IND	15-950-IND
Pink	15-906-IND	15-916-IND	15-960-IND
Orange	15-907-IND	15-917-IND	15-970-IND
Black	15-909-IND	15-919-IND	15-990-IND
Light Blue	15-922-IND	15-932-IND	15-942-IND

BATRIK Superior Face Protection SAVEFACE™ Shields use superior material (7.5 mil thickness PET) to ensure excellent visibility and extended usage due to the design and material quality.

Features

- Flat viewing area that reduces distortion and glare
- Self-shaping foam head piece molds comfortably to forehead
- High quality 7.5 mil thickness PET
- Shields can be used for longer periods due to comfort and improved materials
- Superior anti-fog properties on both sides of the shield
- Wide elastic band 1" (25mm) provides additional comfort
- Ideal safety addition to all pandemic kits
- Latex free
- Level 3 protection drape

REDUCE YOUR SHIELD USAGE DUE TO OUR IMPROVED DESIGN, QUALITY AND COMFORT!

SaveFace™ Shield



Shield Length
9" / 23 cm

Item # SF-100

Description

SAVEFACE™ Shield Length – 9" (23 cm)

Packaging**

100 / Box

SafeFace™ Shield with Drape



Shield Length
9" / 23 cm

Item # SF-180

Description

SAVEFACE™ Total Length 15" (38 cm)

Packaging

100 / Box

Drape Length
6" / 15 cm

*Fluid drape barrier includes closing system below the chin for even greater protection
** More units per box reduces transport costs and waste.*

ver 06/04/2019



BATRIK Scrub Pads are ideal for pre-cleaning your surgical instruments without leaving any metal residue that often occurs when using steel wool or metal sponges. The abrasiveness is perfectly balanced to provide you with efficient cleaning without damaging your instruments. Its fine fibers allow the users to **polish, clean** and **debur** all type of material surfaces.

- Abrasive Material: Aluminum Oxide
- Grit: 400
- Style: Flexible
- Pad material: Non-Woven
- Color: Maroon
- Grade: Very Fine
- Type: Hand Pad



Scrub Pads


Item #	Description	Size	Packaging
60-46050	Scrub Pads for Instruments	4" x 6" (10 cm x 15 cm)	150 units / Box (50 bags of 3 units each per box)
60-46003	Scrub Pads for Instruments	4" x 6" (10 cm x 15 cm)	3 u. / bag

ver 06/20/2019

NEW!

BATRIK Enzymatic Detergents - PrAmLi Power³ Enzymes™

The **PrAmLi Power³ Enzymes™** is a unique blend of three powerful enzymes, along with a surfactant to provide the ultimate in cleaning efficacy.

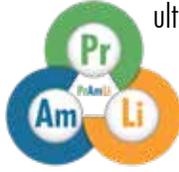

All **PrAmLi Power³ Enzymes™** detergents are multi-enzyme pre-soak cleaners for medical equipment and instrumentation. Our proprietary unique blend of enzymes is effective in the removal of all proteins, blood, fatty tissue and more. It will also remove stains in surgical scrubs, linens and medical clothing.

Products offered with PrAmLi blend:

- Detergents
- Concentrate
- Ready to Use Foaming Spray (**Your First**™) **see p. 44**
- Enzymatic Sponges (**Your First**™) **see p. 46**
- Lubricants (**PrAmLi Power³ Lube™**) **see p. 43**

LOW FOAM DETERGENT

Item #	Packaging	Dilution
		Detergent / Water
SS-21-1001	1 x 3.8 L (1 gallon)	4 to 5 mL / L (0.5 to 1 oz / gallon)
SS-21-1004	4 x 3.8 L (4 x 1 gallon)	4 to 5 mL / L (0.5 to 1 oz / gallon)
SS-21-1025	2 x 9.5 L (2 x 2.5 gallon)	4 to 5 mL / L (0.5 to 1 oz / gallon)



MILD FOAM DETERGENT

Item #	Packaging	Dilution
		Detergent / Water
SS-21-2001	1 x 3.8 L (1 gallon)	4 to 5 mL / L (0.5 to 1 oz / gallon)
SS-21-2004	4 x 3.8 L (4 x 1 gallon)	4 to 5 mL / L (0.5 to 1 oz / gallon)



SUPER CONCENTRATE

Item #	Packaging	Dilution
		Detergent / Water
SS-21-3001	1 x 3.8 L (1 gallon)	1 to 2 mL / L (0.125 to 0.25 oz / gallon)
SS-21-3004	4 x 3.8 L (4 x 1 gallon)	1 to 2 mL / L (0.125 to 0.25 oz / gallon)
SS-21-3025	2 x 9.5 L (2 x 2.5 gallon)	1 to 2 mL / L (0.125 to 0.25 oz / gallon)



ver 01/28/2020



BATRIK Lubricants - PrAmLi Power³ LUBE™

NEW!



PrAmLi Power³ LUBE™ is uniquely formulated which contains an all natural corrosion inhibitor. Specially formulated to lubricate surgical instruments, hinges, box locks and general equipment during autoclaving.

- Forms a protective coating which can extend the instruments useful life
- Bacteriostatic
- Neutral pH 6.5 + 7.5
- Prevents corrosion with its anticorrosive agent

INSTRUMENT LUBRICANT

Item #	Packaging	Dilution Lubricant / Water
SS-21-6001	1 gallon (1 x 3.8 L)	1 to 4 mL / L (0.125 to 0.5 oz / gallon)
SS-21-6004	4 x 1 gallon (4 x 3.8 L) 2 x 2.5 gallon	1 to 4 mL / L (0.125 to 0.5 oz / gallon)
SS-21-6025	(2 x 9.5 L)	1 to 4 mL / L (0.125 to 0.5 oz / gallon)



LUBRICANT - READY TO USE SPRAY

Item #	Packaging
SS-21-6032	1 x 946 mL spray bottle (1 x 32 oz)
SS-21-6012	12 x 946 mL spray bottle (12 x 32 oz)



ver 04/14/2020

BATRIK Foaming Spray is the ideal **Your First Check** to the instrument pre-cleaning process. The thick foam spray breaks down all bioburden to allow for effective cleaning and sterilization. Our Enzyme Foaming Spray is **Your First ™**.



Also with **PrAmLi Power³ Enzymes™** unique blend



ver 04/14/2020

Ready to Use Foaming Spray

Item #	Packaging
SS-21-4001	1 x 946 mL spray bottle (1 x 32 oz)
SS-21-4012	12 x 946 mL spray bottle (12 x 32 oz)

BATRIK Your First ™ Sponges are specifically fabricated for use with all medical instrumentation. Unlike other colored sponges, our sponges are white to better evidence the effectiveness of the cleaning process. Dark colors easily hide the contaminants, where with a white sponge you can visually evaluate the pre-clean efficacy.

Available in two formats, Flat and Tubular:

- 100% medical grade high-density foam
- Will not flake or scratch expensive equipment
- High absorbency as it retains more solution than conventional sponges or paper wipes
- Compatible with enzymatic or disinfectant solutions
- Eco-friendly

See our line of enzymatic sponges!



Can also be used as a bumper to protect your scopes!

Non-Enzymatic Sponges

Item #	Description	Packaging
17-100	Tubular Sponge / White Dispenser Box	25 Units / Box – 4 Boxes / Case
17-101	Tubular / Bulk	560 Units / Case
17-200	Flat / White Dispenser Box	50 Units / Box – 4 Boxes / Case
17-201	Flat / Bulk	1,000 Units / Case

Note: The white sponge may at times appear slightly yellowed. This is only due to UV exposure, not a reflection of its age.

ver 04/14/2020

BATRIK Bedside Sponge Kits With Non-Enzymatic Detergent are **Your First ™** in the medical instrument cleaning process.



Bed Side Kits with non-enzymatic detergent

Item #	Description	Packaging
717-122	Dry Tubular Sponge in large freestanding pouch to use as a basin and simply add 500 mL of water.	15 units / box, 4 boxes / 60 case
717-222	Dry Flat Sponge in large freestanding pouch to use as a basin and simply add 500 mL of water.	25 units / box, 4 boxes / 100 case
717-322	Dry Wave Sponge in large freestanding pouch to use as a basin and simply add 500 mL of water.	Coming Soon
717-133	Tubular Sponge bathing in 500 mL of detergent solution. Ready to Use.	24 units / case
717-233	Flat Sponge bathing in 500 mL of detergent solution. Ready to Use.	24 units / case
717-333	Wave Sponge bathing in 500 mL of detergent solution. Ready to Use.	Coming Soon
717-144	Dual Purpose Tubular Sponge in large freestanding pouch soaked in 100 mL of detergent solution. Ready to Use.	12 units / box, 4 boxes / 48 case
717-244	Dual Purpose Flat Sponge in large freestanding pouch soaked in 100 mL of detergent solution. Ready to Use.	12 units / box, 4 boxes / 48 case
717-344	Dual Purpose Wave Sponge in large freestanding pouch soaked in 100 mL of detergent solution. Ready to Use.	Coming Soon
717-155	Tubular Sponge bathing in 220 mL of detergent solution. Ready to Use.	24 units / case
717-255	Flat Sponge bathing in 220 mL of detergent solution. Ready to Use.	24 units / case
717-355	Wave Sponge bathing in 220 mL of detergent solution. Ready to Use.	Coming Soon

ver-04/14/2020

BATRIK Bedside Enzymatic Sponge Kits are **Your First ™** in the endoscope cleaning process. Our sponges are all pre-treated with **PrAmLi Power³ Enzymes™** to provide the most effective pre-clean or surface clean solution. This unique three-enzyme blend provides you with confident and effective cleaning. The tubular sponge is ideal for all endoscopes.



Enzymatic Bedside Kits

Item #	Description	Packaging
17-111	Dry Tubular Sponge. Moisten in any 250-500 mL basin to activate the enzymatic concentrate.	25 units / box, 4 boxes / 100 case
17-211	Dry Flat Sponge. Moisten in any 250-500 mL basin to activate the enzymatic concentrate.	50 units / box, 4 boxes / 200 case
17-311	Dry Wave Sponge. Moisten in any 250-500 mL basin to activate the enzymatic concentrate.	Coming Soon
17-122	Dry Tubular Sponge in large freestanding pouch to use as a basin and simply add 500 mL of water.	15 units / box, 4 boxes / 60 case
17-222	Dry Flat Sponge in large freestanding pouch to use as a basin and simply add 500 mL of water.	25 units / box, 4 boxes / 100 case
17-322	Dry Wave Sponge in large freestanding pouch to use as a basin and simply add 500 mL of water.	Coming Soon
17-133	Tubular Sponge bathing in 500 mL of enzymatic solution. Ready to Use.	24 units / case
17-233	Flat Sponge bathing in 500 mL of enzymatic solution. Ready to Use.	24 units / case
17-333	Wave Sponge bathing in 500 mL of enzymatic solution. Ready to Use.	Coming Soon
17-144	Dual Purpose Tubular Sponge in large freestanding pouch soaked in 100 mL of enzymatic solution. Ready to Use.	12 units / box, 4 boxes / 48 case
17-244	Dual Purpose Flat Sponge in large freestanding pouch soaked in 100 mL of enzymatic solution. Ready to Use.	12 units / box, 4 boxes / 48 case
17-344	Dual Purpose Wave Sponge in large freestanding pouch soaked in 100 mL of enzymatic solution. Ready to Use.	Coming Soon
17-155	Tubular Sponge bathing in 220 mL of enzymatic solution. Ready to Use.	24 units / case
17-255	Flat Sponge bathing in 220 mL of enzymatic solution. Ready to Use.	24 units / case
17-355	Wave Sponge bathing in 220 mL of enzymatic solution. Ready to Use.	Coming Soon



BATRIK FIRST Wipes are Cleanroom wipes destined to do the first wipe-down of an endoscope after a procedure. They are made from hydroentangled non-woven material composed of 45% polyester and 55% cellulose. These wipes are manufactured and packaged in a cleanroom environment.

Features

- Compatible with virtually any solvent or solution for cleaning and wiping (such as **PrAmLi Power³ Enzymes™**)
- Double-bag packaging
- Production process guarantees extremely low particle generation and chemical extractables

Cleanroom Environment

- Can be used for a cleanroom environment of ISO Class 5 through 8



FIRST Wipes

Item #	Description	Size	Packaging
17-A0909	Nonwoven Cleanroom Wipes	9" x 9" (23 cm x 23 cm)	300 u. / bag 14 bags / case
17-A1212	Nonwoven Cleanroom Wipes	12" x 12" (31 cm x 31 cm)	150 u. / bag 10 bags / case

ver 04/14/2020



BATRIK Transport Bags for Endoscope and Medical Instruments are thoughtfully designed to eliminate touch-points during endoscope transport from point-of-care to the reprocessing area, and back, focused on safety of patients and medical staff.

Two separate disposable clear plastic bags were designed for that purpose:

- Made of recyclable LDPE material
- Excellent resistance to dilute and concentrated acids, alcohols, bases and esters*
- Can withstand high temperatures of up to 90 °C (194 °F)
- Is characterized by being quite flexible, highly resilient, and tough
- Good resistance to aldehydes, ketones and vegetable-based oils*
- Latex and DEHP-free
- Tamper-evident seal for added safety, also aids in keeping the instrument moist to avoid drying of bio-debris
- Clearly marked bags for immediate identification of clean and used scopes with the help of color-coded visuals
- Users can write important information directly on the bag in spaces marked with visual icons
- Made in Canada

Transport Bags

Item #	Description	Size Inches / cm	Packaging
17-HAP24	For Clean Instruments	24" x 22" / 61.0 cm x 56.0 cm	200 Units / Case
17-SAD24	For Soiled Instruments	24" x 22" / 61.0 cm x 56.0 cm	200 Units / Case
17-HAP30	For Clean Instruments	30.5" x 26" / 77.0 cm x 66.0 cm	200 Units / Case
17-SAD30	For Soiled Instruments	30.5" x 26" / 77.0 cm x 66.0 cm	200 Units / Case

For Soiled instruments and accessories that require reprocessing



Item #17-SAD24

Length
24" / 61.0 cm

Item #17-SAD30

Length
30.5" / 77.0 cm

Item #17-SAD24
Width 22" / 56.0 cm

Item #17-SAD30
Width 26" / 66.0 cm

* Reference on file.

For clean instruments and accessories ready for use



Item #17-HAP24

Length
24" / 61.0 cm

Item #17-HAP30

Length
30.5" / 77.0 cm

Item #17-HAP24
Width 22" / 56.0 cm

Item #17-HAP30
Width 26" / 66.0 cm

Using **Your First ™ SEALABLE** transport bags reduces user error, prevents cross-contamination and promotes patient and medical staff safety.

BATRIK Single Use and Reusable Wire Style Brushes - ARC™

Brush Selection Tips:

PART 1: Instrument Specifications

A. What is the maximum length of the instrument(s) to be cleaned?

Add at least 2" (5cm) to that length, so you have an adequate length to work with.

B. What is the inner diameter of the instrument?

Add at most 0.125" (3mm) to determine the diameter of the bristle to ensure effective cleaning. The bristle diameter has to slightly exceed the inner tube diameter. Bristles that are too large will bend too much and not clean effectively.

PART 2: Brush Selection

Now that you have the basic information for your required brush, simply find the section in our catalogue that best corresponds to the overall length required as determined in PART 1. Now select the brush that best corresponds to the diameter determined in PART 2.

As a rule, if your present brush is 8" (20.32cm) in length and we only offer 10" (25.4cm) length with your required diameter, that is entirely acceptable.

In many cases, the BATRIK Brush will have more bristles on the brush. As an example, many 5mm (0.19") diameter brushes are only 1" (2.5cm) in length. At BATRIK, we have strived to provide the greatest amount of bristles to ensure the most effective cleaning.

At Last... Customizable Wire Brushes!

A unique aspect of single-use ARC™ Wire Brushes is that they can be customized to any length:

- Up to 2000mm for Wire Brushes on plastic semi-rigid sheaths, rigid brushes from 5 to 20mm in diameter
- No limit in length for Wire Brushes on flexible plastic sheaths, flexible brushes from 2-5mm in diameter

BATRIK Arc Brushes

ARC™ Wire Brushes are guaranteed to last longer than twisted wire and preserve the integrity of your valuable instruments.



FLEXIBLE
BRUSH & SHEATH

RIGID
BRUSH & SHEATH

12" / 30cm Single-End Wire Brush

	Single Use Item #	Steamable Item #	Color	Total Length inches / cm	Brush (Ø) inches / mm	Brush Length inches / mm	Sheath (Ø) inches / mm
FLEXIBLE BRUSH & SHEATH	41-31502-S	41-31502	White	12.401" / 31.5cm	0.078" / 2mm	0.59" / 15mm	0.39" / 1mm
	41-32503-S	41-32503	Yellow	12.795" / 32.5cm	0.118" / 3mm	0.98" / 25mm	0.39" / 1mm
	411-32504-S	411-32504	Violet	12.795" / 32.5cm	0.157" / 4mm	0.98" / 25mm	0.61" / 1.55mm
	411-32505-S	411-32505	Grey	12.795" / 32.5cm	0.19" / 5mm	0.98" / 25mm	0.68" / 1.75mm
RIGID BRUSH & SHEATH	41-35005-S	41-35005	White	13.779" / 35cm	0.19" / 5mm	1.968" / 50mm	0.118" / 3mm
	41-35006-S	41-35006	Green	13.779" / 35cm	0.23" / 6mm	1.968" / 50mm	0.118" / 3mm
	41-35008-S	41-35008	Blue	13.779" / 35cm	0.31" / 8mm	1.968" / 50mm	0.157" / 4mm
	41-35010-S	41-35010	Red	13.779" / 35cm	0.39" / 10mm	1.968" / 50mm	0.157" / 4mm
	41-35012-S	41-35012	Yellow	13.779" / 35cm	0.49" / 12.5mm	1.968" / 50mm	0.157" / 4mm
	41-35015-S	41-35015	Green	13.779" / 35cm	0.59" / 15mm	1.968" / 50mm	0.228" / 5.8mm
	41-35020-S	41-35020	Blue	13.779" / 35cm	0.78" / 20mm	1.968" / 50mm	0.228" / 5.8mm

Packaging: Single Use: Bag of 20 / Steamable: Bag of 3



18" / 45cm Single-End Wire Brush

	Single Use Item #	Steamable Item #	Color	Total Length Inches / cm	Brush (Ø) Inches / mm	Brush Length Inches / mm	Sheath (Ø) Inches / mm
FLEXIBLE BRUSH & SHEATH	41-46502-S	41-46502	White	18.307" / 46.5cm	0.078" / 2mm	0.59" / 15mm	0.39" / 1mm
	41-47503-S	41-47503	Yellow	18.7" / 47.5cm	0.118" / 3mm	0.98" / 25mm	0.39" / 1mm
	411-47504-S	411-47504	Violet	18.7" / 47.5cm	0.157" / 4mm	0.98" / 25mm	0.61" / 1.55mm
	411-47505-S	411-47505	Grey	18.7" / 47.5cm	0.19" / 5mm	0.98" / 25mm	0.68" / 1.75mm
RIGID BRUSH & SHEATH	41-50005-S	41-50005	White	19.685" / 50cm	0.19" / 5mm	1.968" / 50mm	0.118" / 3mm
	41-50006-S	41-50006	Green	19.685" / 50cm	0.23" / 6mm	1.968" / 50mm	0.118" / 3mm
	41-50008-S	41-50008	Blue	19.685" / 50cm	0.31" / 8mm	1.968" / 50mm	0.157" / 4mm
	41-50010-S	41-50010	Red	19.685" / 50cm	0.39" / 10mm	1.968" / 50mm	0.157" / 4mm
	41-50012-S	41-50012	Yellow	19.685" / 50cm	0.49" / 12.5mm	1.968" / 50mm	0.157" / 4mm
	41-50015-S	41-50015	Green	19.685" / 50cm	0.59" / 15mm	1.968" / 50mm	0.228" / 5.8mm
	41-50020-S	41-50020	Blue	19.685" / 50cm	0.78" / 20mm	1.968" / 50mm	0.228" / 5.8mm

Packaging: Single Use: Bag of 20 / Steamable: Bag of 3

24" / 60cm Single-End Wire Brush

	Single Use Item #	Steamable Item #	Color	Total Length Inches / cm	Brush (Ø) Inches / mm	Brush Length Inches / mm	Sheath (Ø) Inches / mm
FLEXIBLE BRUSH & SHEATH	41-61502-S	41-61502	White	24.212" / 61.5cm	0.078" / 2mm	0.59" / 15mm	0.39" / 1mm
	41-62503-S	41-62503	Yellow	24.6" / 62.5cm	0.118" / 3mm	0.98" / 25mm	0.39" / 1mm
	411-62504-S	411-62504	Violet	24.6" / 62.5cm	0.157" / 4mm	0.98" / 25mm	0.61" / 1.55mm
	411-62505-S	411-62505	Grey	24.6" / 62.5cm	0.19" / 5mm	0.98" / 25mm	0.68" / 1.75mm
RIGID BRUSH & SHEATH	41-65005-S	41-65005	White	25.59" / 65cm	0.19" / 5mm	1.968" / 50mm	0.118" / 3mm
	41-65006-S	41-65006	Green	25.59" / 65cm	0.23" / 6mm	1.968" / 50mm	0.118" / 3mm
	41-65008-S	41-65008	Blue	25.59" / 65cm	0.31" / 8mm	1.968" / 50mm	0.157" / 4mm
	41-65010-S	41-65010	Red	25.59" / 65cm	0.39" / 10mm	1.968" / 50mm	0.157" / 4mm
	41-65012-S	41-65012	Yellow	25.59" / 65cm	0.49" / 12.5mm	1.968" / 50mm	0.157" / 4mm
	41-65015-S	41-65015	Green	25.59" / 65cm	0.59" / 15mm	1.968" / 50mm	0.228" / 5.8mm
	41-65020-S	41-65020	Blue	25.59" / 65cm	0.78" / 20mm	1.968" / 50mm	0.228" / 5.8mm


Packaging: Single Use: Bag of 20 / Steamable: Bag of 3

30" / 75cm Single-End Wire Brush


	Single Use Item #	Steamable Item #	Color	Total Length Inches / cm	Brush (Ø) Inches / mm	Brush Length Inches / mm	Sheath (Ø) Inches / mm
FLEXIBLE BRUSH & SHEATH	41-76502-S	41-76502	White	30.118" / 76.5cm	0.078" / 2mm	0.59" / 15mm	0.39" / 1mm
	41-77503-S	41-77503	Yellow	30.511" / 77.5cm	0.118" / 3mm	0.98" / 25mm	0.39" / 1mm
	411-77504-S	411-77504	Violet	30.511" / 77.5cm	0.157" / 4mm	0.98" / 25mm	0.61" / 1.55mm
	411-77505-S	411-77505	Grey	30.511" / 77.5cm	0.19" / 5mm	0.98" / 25mm	0.68" / 1.75mm
RIGID BRUSH & SHEATH	41-80005-S	41-80005	White	31.496" / 80cm	0.19" / 5mm	1.968" / 50mm	0.118" / 3mm
	41-80006-S	41-80006	Green	31.496" / 80cm	0.23" / 6mm	1.968" / 50mm	0.118" / 3mm
	41-80008-S	41-80008	Blue	31.496" / 80cm	0.31" / 8mm	1.968" / 50mm	0.157" / 4mm
	41-80010-S	41-80010	Red	31.496" / 80cm	0.39" / 10mm	1.968" / 50mm	0.157" / 4mm
	41-80012-S	41-80012	Yellow	31.496" / 80cm	0.49" / 12.5mm	1.968" / 50mm	0.157" / 4mm
	41-80015-S	41-80015	Green	31.496" / 80cm	0.59" / 15mm	1.968" / 50mm	0.228" / 5.8mm
	41-80020-S	41-80020	Blue	31.496" / 80cm	0.78" / 20mm	1.968" / 50mm	0.228" / 5.8mm

Packaging: Single Use: Bag of 20 / Steamable: Bag of 3

FLEXIBLE

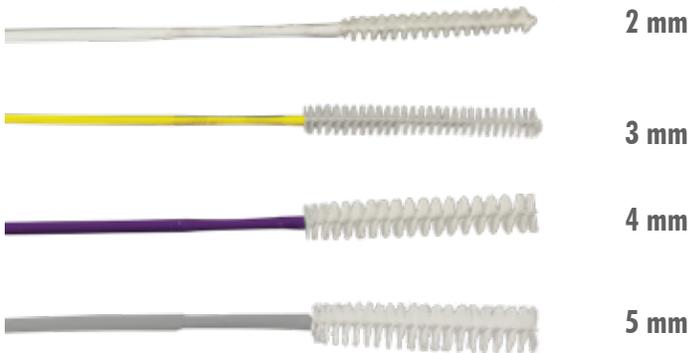

RIGID

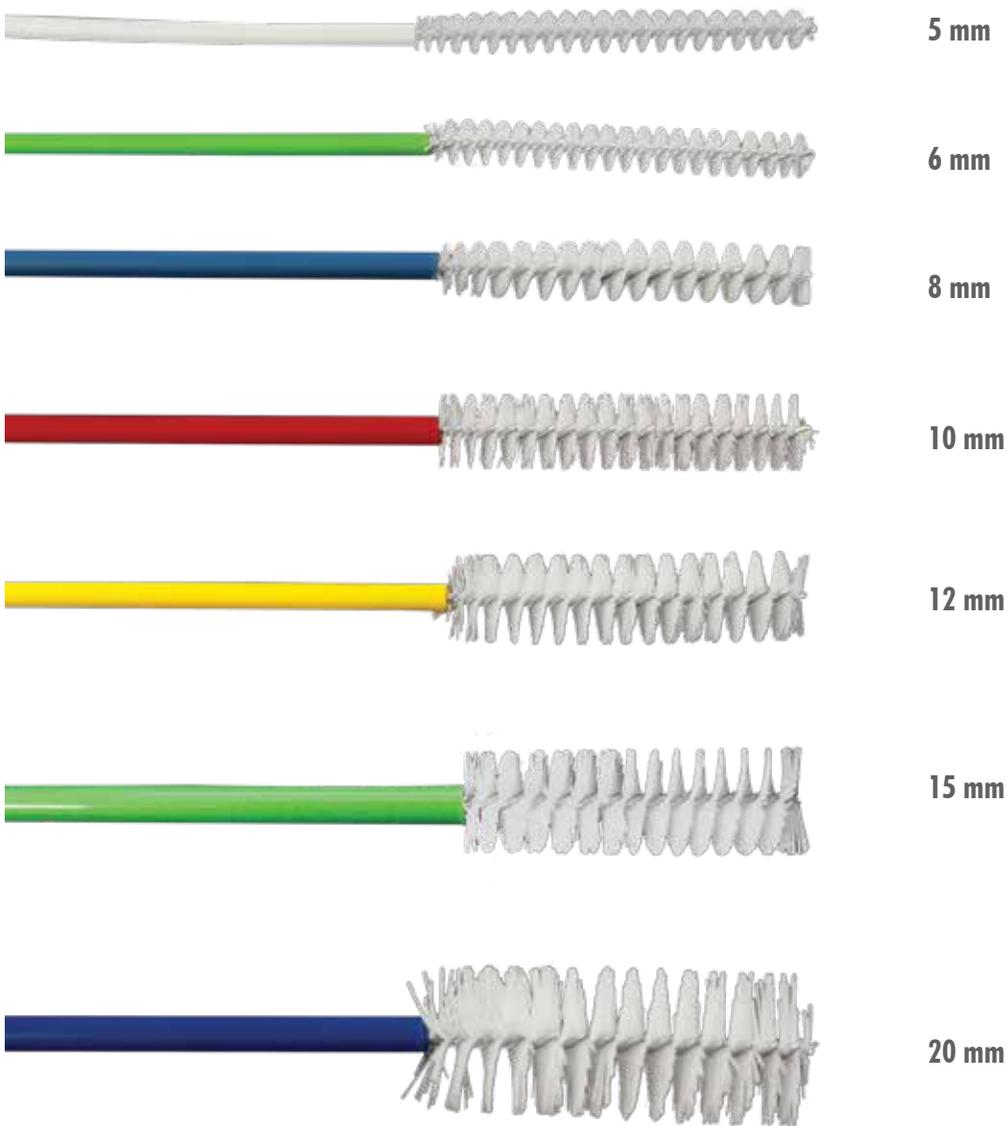

vet04/14/2020


BATRIK Reusable Orthopedic Rasps Cleaning Brush

NEW!
Antibacterial

This reusable brush is used to clean orthopedic rasps and other similar instruments in orthopaedic surgery.

Features

- Antibacterial polyamide bristles
- Ergonomic steamable thermoplastic handle
- Reusable
- Long lasting and can withstand many sterilization cycles
- Latex free

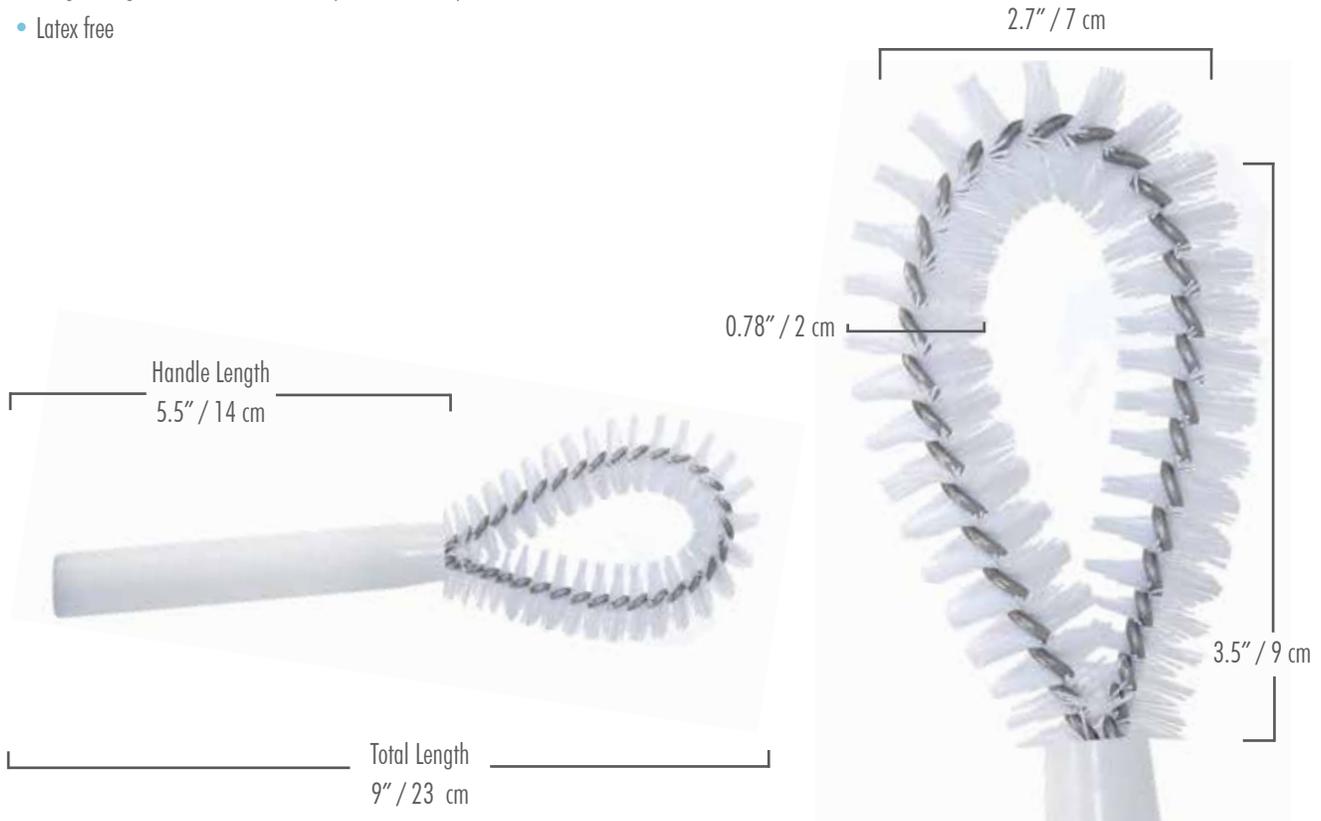

Orthopedic Rasps Cleaning Brush

Item #	Total Length Inches / cm	Handle Length Inches / cm	Brush Size Inches / cm	Packaging
50-RASP	9" / 23 cm	5.5" / 14 cm	0.78" x 2.7" x 3.5" 2 cm x 7 cm x 9 cm	3 Units / Bag

ANTIBACTERIAL

The bristles are made of polyamide fibers which prevent or slow bacterial growth between two cleanings or brush sterilizations.

ver 04/14/2020

NEW!
1.5 mm Brush

BATRIK Reusable Stainless Steel ARC™ Micro Brushes

Our Unique 1.5 mm Stainless Steel Brush is Now Available!

One of the smallest brushes in the industry: They are designed to clean your channels with an inner diameter ranging from 0.5 mm to 1.40 mm.

It is highly recommended using the **pull-through technique** (insert the blunt, non-brush end first).

Features

- Flexible stainless steel sheath - will not bend like the regular metal twisted wire brushes
- Polyamide bristles
- Reusable*
- Packaging: 6 units / Tube



1.5 mm Brush

Item #	Total Length Inches / cm	Brush Size (Ø) Inches / mm	Brush Length Inches / mm	Sheath (Ø) Inches / mm
35-03015	12" / 30 cm	0.06" / 1.5 mm	0.6" / 15 mm	0.02" / 0.40 mm
35-04515	18" / 45 cm	0.06" / 1.5 mm	0.6" / 15 mm	0.02" / 0.40 mm

*** Due to the delicate nature of these brushes, although all the material used in the fabrication can be re-processed repeatedly (as with any other wire brush), any re-use is conditional as to the handling and application of this brush.**

ver 04/14/2020

BATRIK Reusable Asymmetric Double ARC™ Wire Brushes**14" / 35cm Asymmetric Double-End Wire Brush**

Item #	Color	Length Inches / cm	Brush 1 (Ø) Inches / mm	Brush 2 (Ø) Inches / mm	Brush Length Inches / mm	Sheath (Ø) Inches / mm
40-35056	White	13.779" / 35cm	0.19" / 5 mm	0.23" / 6mm	1.968" / 50mm	0.118" / 3mm
40-35810	Blue	13.779" / 35cm	0.31" / 8mm	0.39" / 10mm	1.968" / 50mm	0.157" / 4mm
40-35112	Yellow	13.779" / 35cm	0.39" / 10mm	0.49" / 12.5mm	1.968" / 50mm	0.157" / 4mm
40-35152	Green	13.779" / 35cm	0.59" / 15mm	0.78" / 20mm	1.968" / 50mm	0.228" / 5.8mm

Packaging: Bag of 3**20" / 50cm Asymmetric Double-End Wire Brush**


Item #	Color	Length Inches / cm	Brush 1 (Ø) Inches / mm	Brush 2 (Ø) Inches / mm	Brush Length Inches / mm	Sheath (Ø) Inches / mm
40-50056	White	19.685" / 50cm	0.19" / 5 mm	0.23" / 6mm	1.968" / 50mm	0.118" / 3mm
40-50810	Blue	19.685" / 50cm	0.31" / 8mm	0.39" / 10mm	1.968" / 50mm	0.157" / 4mm
40-50112	Yellow	19.685" / 50cm	0.39" / 10mm	0.49" / 12.5mm	1.968" / 50mm	0.157" / 4mm
40-50152	Green	19.685" / 50cm	0.59" / 15mm	0.78" / 20mm	1.968" / 50mm	0.228" / 5.8mm

Packaging: Bag of 3**26" / 65cm Asymmetric Double-End Wire Brush**

Item #	Color	Length Inches / cm	Brush 1 (Ø) Inches / mm	Brush 2 (Ø) Inches / mm	Brush Length Inches / mm	Sheath (Ø) Inches / mm
40-65056	White	25.59" / 65cm	0.19" / 5 mm	0.23" / 6mm	1.968" / 50mm	0.118" / 3mm
40-65810	Blue	25.59" / 65cm	0.31" / 8mm	0.39" / 10mm	1.968" / 50mm	0.157" / 4mm
40-65112	Yellow	25.59" / 65cm	0.39" / 10mm	0.49" / 12.5mm	1.968" / 50mm	0.157" / 4mm
40-65152	Green	25.59" / 65cm	0.59" / 15mm	0.78" / 20mm	1.968" / 50mm	0.228" / 5.8mm

Packaging: Bag of 3**31" / 80cm Asymmetric Double-End Wire Brush**

Item #	Color	Length Inches / cm	Brush 1 (Ø) Inches / mm	Brush 2 (Ø) Inches / mm	Brush Length Inches / mm	Sheath (Ø) Inches / mm
40-80056	White	31.496" / 80cm	0.19" / 5 mm	0.23" / 6mm	1.968" / 50mm	0.118" / 3mm
40-80810	Blue	31.496" / 80cm	0.31" / 8mm	0.39" / 10mm	1.968" / 50mm	0.157" / 4mm
40-80112	Yellow	31.496" / 80cm	0.39" / 10mm	0.49" / 12.5mm	1.968" / 50mm	0.157" / 4mm
40-80152	Green	31.496" / 80cm	0.59" / 15mm	0.78" / 20mm	1.968" / 50mm	0.228" / 5.8mm

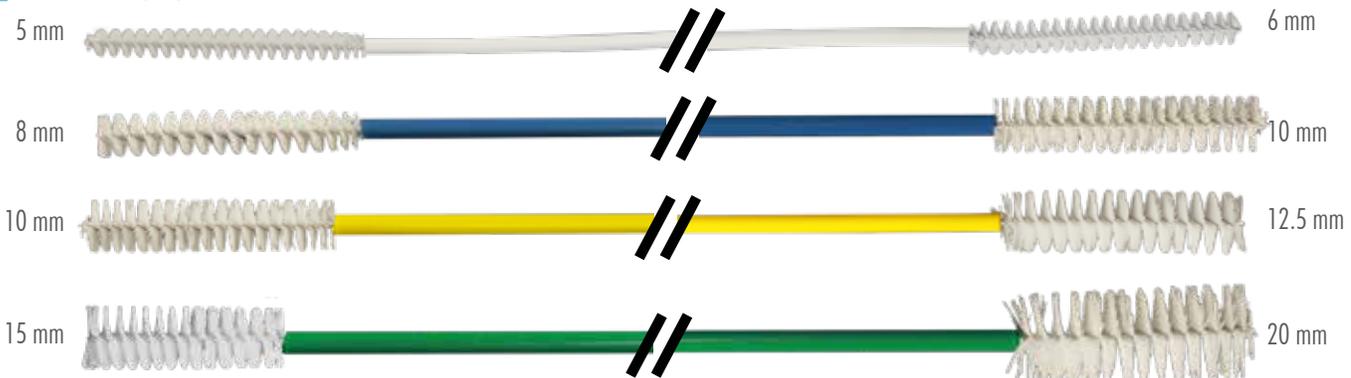Packaging: Bag of 3

ver 05/05/2020

BATRIK Reusable Symmetric Double ARC™ Wire Brushes

14" / 35cm Symmetric Double-End Wire Brushes

Item #	Color	Length Inches / cm	Brush 1 (Ø) Inches / mm	Brush Length Inches / mm	Sheath (Ø) Inches / mm
40-35005	White	13.779" / 35cm	0.19" / 5 mm	1.968" / 50mm	0.118" / 3mm
40-35006	Green	13.779" / 35cm	0.23" / 6mm	1.968" / 50mm	0.118" / 3mm
40-35008	Blue	13.779" / 35cm	0.31" / 8mm	1.968" / 50mm	0.157" / 4mm
40-35010	Red	13.779" / 35cm	0.39" / 10mm	1.968" / 50mm	0.157" / 4mm
40-35012	Yellow	13.779" / 35cm	0.49" / 12.5mm	1.968" / 50mm	0.157" / 4mm
40-35015	Green	13.779" / 35cm	0.59" / 15mm	1.968" / 50mm	0.228" / 5.8mm
40-35020	Blue	13.779" / 35cm	0.78" / 20mm	1.968" / 50mm	0.228" / 5.8mm

Packaging: Bag of 3

20" / 50cm Symmetric Double-End Wire Brushes

Item #	Color	Length Inches / cm	Brush 1 (Ø) Inches / mm	Brush Length Inches / mm	Sheath (Ø) Inches / mm
40-50005	White	19.685 / 50cm	0.19" / 5 mm	1.968" / 50mm	0.118" / 3mm
40-50006	Green	19.685 / 50cm	0.23" / 6mm	1.968" / 50mm	0.118" / 3mm
40-50008	Blue	19.685 / 50cm	0.31" / 8mm	1.968" / 50mm	0.157" / 4mm
40-50010	Red	19.685 / 50cm	0.39" / 10mm	1.968" / 50mm	0.157" / 4mm
40-50012	Yellow	19.685 / 50cm	0.49" / 12.5mm	1.968" / 50mm	0.157" / 4mm
40-50015	Green	19.685 / 50cm	0.59" / 15mm	1.968" / 50mm	0.228" / 5.8mm
40-50020	Blue	19.685 / 50cm	0.78" / 20mm	1.968" / 50mm	0.228" / 5.8mm

Packaging: Bag of 3

26" / 65cm Symmetric Double-End Wire Brushes

Item #	Color	Length Inches / cm	Brush 1 (Ø) Inches / mm	Brush Length Inches / mm	Sheath (Ø) Inches / mm
40-65005	White	25.59" / 65cm	0.19" / 5 mm	1.968" / 50mm	0.118" / 3mm
40-65006	Green	25.59" / 65cm	0.23" / 6mm	1.968" / 50mm	0.118" / 3mm
40-65008	Blue	25.59" / 65cm	0.31" / 8mm	1.968" / 50mm	0.157" / 4mm
40-65010	Red	25.59" / 65cm	0.39" / 10mm	1.968" / 50mm	0.157" / 4mm
40-65012	Yellow	25.59" / 65cm	0.49" / 12.5mm	1.968" / 50mm	0.157" / 4mm
40-65015	Green	25.59" / 65cm	0.59" / 15mm	1.968" / 50mm	0.228" / 5.8mm
40-65020	Blue	25.59" / 65cm	0.78" / 20mm	1.968" / 50mm	0.228" / 5.8mm

Packaging: Bag of 3

31" / 80cm Symmetric Double-End Wire Brushes

Item #	Color	Length Inches / cm	Brush 1 (Ø) Inches / mm	Brush Length Inches / mm	Sheath (Ø) Inches / mm
40-80005	White	31.496" / 80cm	0.19" / 5 mm	1.968" / 50mm	0.118" / 3mm
40-80006	Green	31.496" / 80cm	0.23" / 6mm	1.968" / 50mm	0.118" / 3mm
40-80008	Blue	31.496" / 80cm	0.31" / 8mm	1.968" / 50mm	0.157" / 4mm
40-80010	Red	31.496" / 80cm	0.39" / 10mm	1.968" / 50mm	0.157" / 4mm
40-80012	Yellow	31.496" / 80cm	0.49" / 12.5mm	1.968" / 50mm	0.157" / 4mm
40-80015	Green	31.496" / 80cm	0.59" / 15mm	1.968" / 50mm	0.228" / 5.8mm
40-80020	Blue	31.496" / 80cm	0.78" / 20mm	1.968" / 50mm	0.228" / 5.8mm

Packaging: Bag of 3

Ver 05/05/2020



BATRIK In-Brush™ ANTIMICROBIAL Twisted Wire Brushes

Twisted Wire Brushes

Designed for use with suction tubes, our twisted wires have a fine / medium diameter yet are rigid enough to ensure excellent brush cleaning capabilities

Nylon Twisted Wire Brushes - with Antimicrobial Bristles

We offer a wide selection of brushes for most of your instrument cleaning applications. Ideal for large bore, uniquely shaped, and endoscopic instruments.

Now with new protective
Blue Tip for added safety



*Please note that our blue-tipped brushes are currently in transition. You may still receive brushes without the blue tip.

10" / 250 mm Frazier / Baron Suction Tubes Brushes

Item #	Size	Brush Length Inches / mm	Brush (Ø) Inches / mm
20-10079	6 FR	1" / 25.4 mm	0.079" / 2 mm
20-10090	7 FR	1" / 25.4 mm	0.090" / 2.3 mm
20-10106	8 FR	1" / 25.4 mm	0.106" / 2.7 mm
20-10118	9 FR	1" / 25.4 mm	0.118" / 3 mm
20-10129	10 FR	1.5" / 38 mm	0.129" / 3.3 mm
20-10146	11 FR	1.5" / 38 mm	0.0146" / 3.7 mm
20-10157	12 FR	1.5" / 38 mm	0.157" / 4 mm

Packaging: Bag of 3

16" / 406 mm Nylon Twisted Wire Brushes

Item #	Brush Length Inches / mm	Brush (Ø) Inches / mm
23-16197	2.5" / 63.5 mm	0.197" / 5 mm
23-16236	2.5" / 63.5 mm	0.236" / 6 mm
23-16276	3" / 76 mm	0.276" / 7 mm
23-16394	3" / 76 mm	0.394" / 10 mm
23-16472	3" / 76 mm	0.472" / 12 mm
23-16591	3" / 76 mm	0.591" / 15 mm
23-16787	3" / 76 mm	0.787" / 20 mm

Packaging: Bag of 3

Broncho Brushes



Item #	Description	Length Inches / mm	Brush Length Inches / mm
26-06175	Small	6" / 152 mm	1.75" / 44.5 mm
26-08325	Large	8" / 203 mm	3.25" / 82.6 mm

Packaging: Bag of 3

Tracheostomy Brushes

Item #	Description	Length Inches / mm	Brush Length Inches / mm	Brush (Ø) Inches / mm
25-08500	Large 6-12.3 ^{3/4"}	8" / 203 mm	3.75" / 95 mm	0.5" / 12.7 mm
25-37550	Small 000-5.1 ^{7/8"}	4" / 93 mm	1.88" / 48 mm	0.5" / 12.7 mm

Packaging: Bag of 3



12" / 305 mm Nylon Twisted Wire Brushes

Item #	Brush Length Inches / mm	Brush (Ø) Inches / mm
22-12157	3" / 76 mm	0.157" / 4 mm
22-12197	3" / 76 mm	0.197" / 5 mm
22-12236	3" / 76 mm	0.236" / 6 mm
22-12393	3" / 76 mm	0.394" / 10 mm

Packaging: Bag of 3

18" / 457 mm Nylon Twisted Wire Brushes


Item #	Brush Length Inches / mm	Brush (Ø) Inches / mm
19-18010	0.47" / 12 mm	0.039" / 1 mm
19-18015	0.5" / 12 mm	0.06" / 1.6 mm

Packaging: Bag of 3

24" / 610 mm Nylon Twisted Wire Brushes

Item #	Brush Length Inches / mm	Brush (Ø) Inches / mm
24-24098	2.5" / 63.5 mm	0.098" / 2.5 mm
24-24197	2.5" / 63.5 mm	0.197" / 5 mm
24-24276	3" / 76 mm	0.276" / 7 mm
24-24394	3" / 76 mm	0.394" / 10 mm
24-24472	3" / 76 mm	0.472" / 12 mm

Packaging: Bag of 3

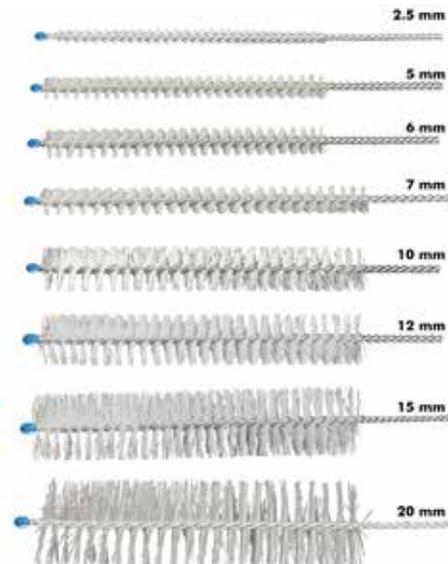

VER 04/14/2020

BATRIK Instrument Cleaning Flat Brushes are specifically designed for hospital use in the pre-clean and preparation of instruments.

- Long, ergonomic designed handle permits ease-of-use
- Extra long length of brush head (60 [± 5] mm in Single-Use and 75 [± 5] mm in Reusable)
- The brushes fibers are transparent, allowing for quick verification of its cleanliness
- The plastic composition used to fabricate the handle will not release any toxic components



BATRIK manufactures four types of flat brushes:

- **EXTRA RIGID** instrument safe polyamide fibers offer steel-like rigidity and allows deep cleaning of soiled instruments, particularly in grooved or uneven surfaces.
- **RIGID** the rigidity of this brush, along with the length of the bristles and handle, improves and simplifies instrument cleaning.
- **FLEXIBLE** This brush is designed for large instruments that require a more gentle cleaning application such as optics, plastics or hard to reach spots which are easier to clean using a flexible bristle.
- **EXTRA FLEXIBLE** The very long bristles eliminate any stress during sanitation, and makes this the ideal brush for fragile instruments (optics, etc.) as well as both pliable and rigid endoscopes. The length of the bristles assist in cleaning areas that would otherwise be difficult to reach with shorter bristles, such as the accessory channel or tip deflection of pliable endoscopes



Extra Rigid Reusable / Autoclavable Flat Brush

Item #	Description	Length – mm	Brush Size – mm
50-40917	Reusable Flat brush (extra rigid)	220mm	75L x 20W x 10H

Packaging: Bag of 3



Rigid Reusable / Autoclavable Flat Brush

Item #	Description	Length – mm	Brush Size – mm
50-40918	Reusable Flat brush (rigid)	220mm	75L x 20W x 15H

Packaging: Bag of 3



Flexible Reusable / Autoclavable Flat Brush

Item #	Description	Length – mm	Brush Size – mm
50-40919	Reusable Flat brush (flexible)	220mm	75L x 20W x 15H

Packaging: Bag of 3



Extra Flexible Reusable / Autoclavable Flat Brush

Item #	Description	Length – mm	Brush Size – mm
50-40920	Reusable Flat brush (extra flexible)	220mm	75L x 20W x 25H

Packaging: Bag of 3



Rigid and Flexible Single Use Flat Brush

Item #	Description	Length – mm	Brush Size – mm
50-40921	Single use flat brush (rigid)	220mm	60L x 20W x 15H
50-40922	Single use flat brush (flexible)	220mm	60L x 20W x 15H

Packaging: Bag of 20

POLYAMIDE FIBERS



Ver 07/25/2018



BATRIK Instrument Cleaning Reusable Precision Brush External cleaning brush is designed for cleaning instrument jaws. It is the ideal brush for the quick removal of difficult debris, the BATRIK Reusable Precision Brush can also easily clean deep into instrument jaws. Although the firm and durable bristles feel like stainless steel, there is absolutely no risk of instrument damage or personal injury.

- Polyamide (TYNEX®) bristles
- Bristles have an extra firm density and superior durability
- Steamable thermoplastic handle
- Ergonomic handle
- Steam / chemical sterilizable
- Latex free



No more
Steel brushes

ELIMINATES THE NEED FOR THE STAINLESS STEEL TOOTHBRUSH-STYLE BRUSH, WHICH CAN DAMAGE VALUABLE INSTRUMENTS!



“**THROW AWAY YOUR WIRE BRUSHES, THIS IS THE BEST OPTION TO EXTEND YOUR INSTRUMENT USEFUL LIFE**”

Reusable Precision Brush

Item #	Description	Packaging
50-40800	Reusable Precision External Cleaning Brush Overall length is 6.29" / 160mm	3 units / Bag

ver 02/14/2017


BATRIK Instrument Cleaning Reusable Jar / Bottle Brush with its long and flexible bristles will maximize the effectiveness of large-surface cleaning. The easy-grip handle also provides excellent maneuverability. Designed specifically for hospital applications.

- High quality polyamide clear bristles
- Easy to identify debris
- NEW ergonomic handle
- Autoclavable
- Latex free




Reusable Jar / Bottle Brush

Item #	Description	Packaging
50-10295	Reusable Jar / Bottle Brush – Autoclavable Overall length is 13.19" / 355mm	3 units / Bag

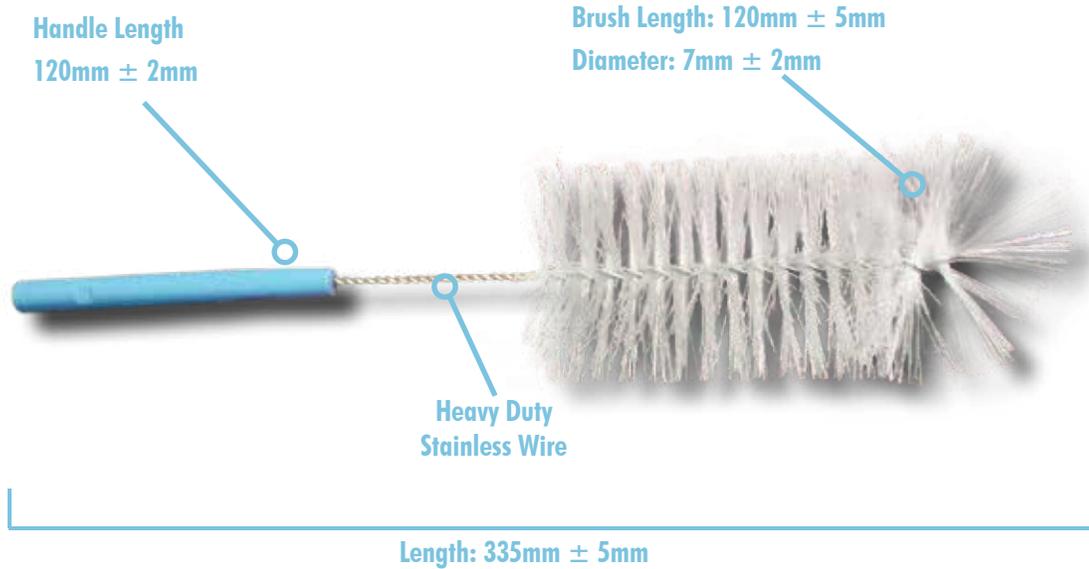

Sterilizer Drain Cleaning Brush

Item #	Description	Packaging
50-10224	Sterilizer Drain Cleaning Brush Overall length is 24" / 609mm	3 units / Bag

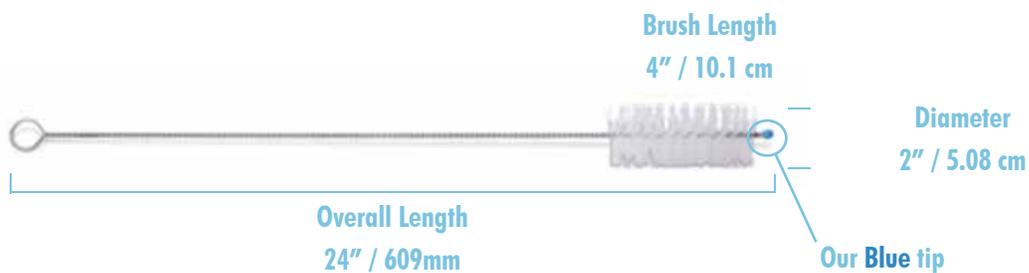

ver 10/30/2018



BATRIK Instrument Cleaning ANTIMICROBIAL Acetabular Reamer Brush - Durable and Reusable Cleaning Brush. This double-ended cleaning brush is manufactured exclusively for cleaning acetabular bone reamers. The "Pig Tail" rounded style circular head brush is designed to effectively clean inside the reamer, while the rigid straight brush on the reverse side is used for cleaning individual cutting holes. Special fusion of the brush heads to the hand piece along with using high quality **antimicrobial** medical grade bristles ensures a durable and reusable effective cleaning tool.

Product Specifications*

- Round brush head diameter: 1.75" (45mm)
- Round brush head length: 2" (50mm)
- Straight brush diameter: 5/16" (8mm) Tapered
- Straight brush length: 1 3/16" (30mm)

*specifications subject to small variances



Brush ±

0.70" x 2" / 18 x 50 mm
Ø 1.75" / 45mm

Handle Length ±

4.72" / 120 mm
Ø 0.47" / 12 mm

Tapered Brush ±

1.18" / 30 mm
Ø 0.3125" / 8 mm

Acetabular Reamer Cleaning Brush

Item #	Description	Packaging
50-ACETAB	Acetabular Reamer Cleaning Brush	3 units / Bag

BATRIK INSTRUMENT CLEANING All Purpose Brushes

Plastic Handle Brushes

We offer a wide selection of brushes for most of your instrument cleaning applications. These brushes are ideal for large bore, uniquely shaped and endoscopic instruments.

In-Brush™ Plastic Handle

Choose either the regular toothbrush style (nylon or stainless) or try our double-ended brush.

*order in volume and save even more

Plastic Single End Nylon – Handle Brush



Item #	Description	Length Inches / mm	Brush Length Inches / mm	Brush (Ø) Inches / mm
28-72550	Single End Nylon	7.25" / 184mm	1.5" / 38mm	0.5" / 12.7mm

Packaging: Bag of 3 - ITEM # 28-72550-50 / Bag of 50

Plastic Single End Stainless Steel – Handle Brush



Item #	Description	Length Inches / mm	Brush Length Inches / mm	Brush (Ø) Inches / mm
29-72550	Single End Stainless Steel	7.25" / 184mm	1.5" / 38mm	0.5" / 12.7mm

Packaging: Bag of 3 - ITEM # 29-72550-50 / Bag of 50

Plastic Single End Brass – Handle Brush



Item #	Description	Length Inches / mm	Brush Length Inches / mm	Brush (Ø) Inches / mm
29-70050	Single End Brass	7.25" / 184mm	1.5" / 38mm	0.5" / 12.7mm

Packaging: Bag of 3 - ITEM # 29-70050-50 / Bag of 50

Plastic Double End Nylon – Handle Brush



Item #	Description	Length Inches / mm	Brush Length Inches / mm	Brush (Ø) Inches / mm
28-69394	Large Brush-end	6.9" / 175mm	1.38" / 35mm	0.39" / 10mm
	Small Brush-end	6.9" / 175mm	1" / 25mm	0.2" / 5mm

Packaging: Bag of 3 - ITEM # 28-69394-50 / Bag of 50

Hand Scrub Nylon Brush – AUTOCLAVABLE

Item #	Length Inches / mm	Brush Length Inches / mm	Brush (Ø) Inches / mm
50-40911	4.09" / 104mm	3.74" / 95mm	1.1" / 28mm



Packaging: Bag of 3 - ITEM # 50-40911-50 / Bag of 50

ver 10/30/2018




BATRIK Instrument Autoclave Cleaning Brush Kits

Regular Autoclave Cleaning Brush Solution Kit includes:

- **1 Hand Brush**
Length: 4.5 x 8" / 11.5 x 20.32cm
- **1 Rounded Corner Brush**
Head size: 5.5 x 10" / 14 x 25.4cm
- **1 Anodized Aluminum Extension Pole**
Size: 2-4ft / 60-122cm

* Brushes are made with FDA approved polyester bristles

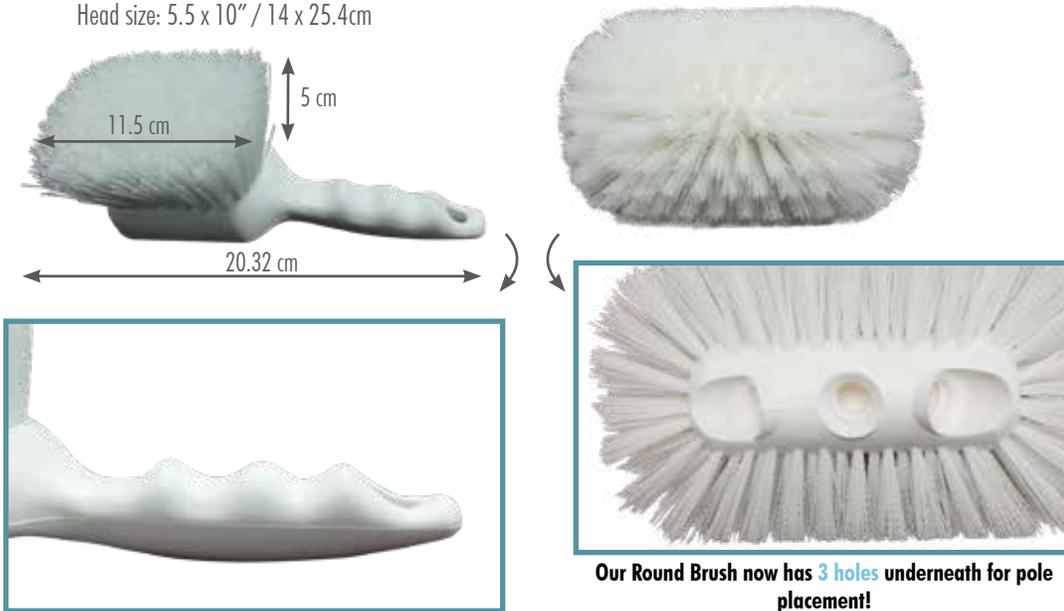

Our new Hand Brush now has an ergonomic handle!

- Better hand comfort
- Longer bristles
- Handle is easier to clean (no cavity)

Our Round Brush now has 3 holes underneath for pole placement!

- **More flexibility** in regards to the brush placement (parallel or perpendicular)
- Hard-to-reach corners will no longer be an issue thanks to the **new brush placement possibilities** and to the **longer brush length** (now 10" instead of 8.75").

*Some brushes may still have only one center hole.



Autoclave Regular Brush Cleaning Kit

Item #	Description
50-CLAVEREG	Autoclave REGULAR Brush Cleaning Kit

Superior Autoclave Cleaning Brush Solution Kit includes:

- **1 Hand Brush**
Length: 4.5 x 8" / 11.5 x 20.32cm
- **1 Rounded Corner Brush**
Head size: 5.5 x 8.75" / 14 x 22.225cm
- **1 Anodized Aluminum Extension Pole**
Size: 2-4ft / 60-122cm
- **2 BATRIK Full Face Shields**

* Brushes are made with FDA approved polyester bristles

Autoclave Superior Brush Cleaning Kit

Item #	Description
50-CLAVESUP	Autoclave SUPERIOR Brush Cleaning Kit

Replacement Components



Hand Brush
50-35900
Packaging: 1



Brush Head
50-55800
Packaging: 1



Extension Pole
50-EXTPOLE
Packaging: 1



BATRIK Face Shields
SF-100
Packaging: 100 / Box

BATRIK Hysterosonography Catheter

Hystero-sonography (also known as sonohysterography) can now be performed safely, simply and atraumatically with the gentle infusion of saline into the uterus with the use of the Batrik H/S Catheter.

Benefits

- Latex Free
- Available in two sizes: 5fr and 7fr
- Minimally invasive and ideal for office use
- Gentle placement sheath included (decreasing the need for a tenaculum) to minimize patient discomfort
- Flexible and easy-to-use
- Total length of 410 mm with a working length of 280 mm



Hysterosonography Catheter

Item no.	Description	Packaging
51-5005	5 FR H/S Catheter	10 / box
51-5007	7 FR H/S Catheter	10 / box

“ WE DEPEND UPON THE RELIABILITY, PRECISION AND FUNCTIONALITY OF THE BATRIK H/S CATHETER IN OUR CLINICAL AND RESEARCH OPERATIONS.


HAVING PERFORMED MORE THAN 100,000 SONOHYSTEROGRAPHY PROCEDURES, WE CAN UNEQUIVOCALLY SAY THAT THIS IS THE BEST HSG CATHETER THAT WE HAVE EVER USED ”



ver 07/28/2016




Improved Imagery and Diagnosis with Sonohysterograms

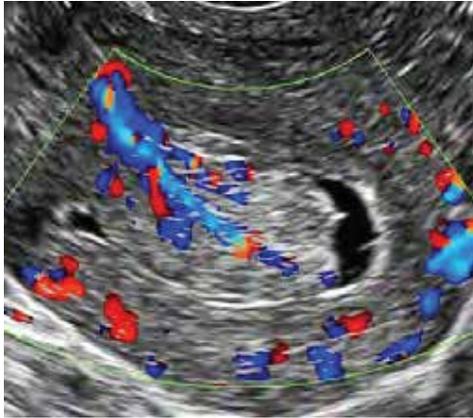


Sonohysterogram showing a polyp with a vascular pedicle in a 28 year old patient with abnormal uterine bleeding and infertility



Sonohysterogram showing a type 0 submucosal broid completely within the uterine cavity

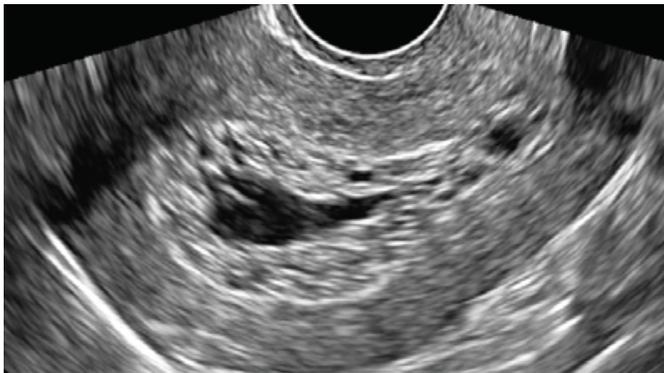

Very clear example of endometrial hyperplasia on sonohysterogram

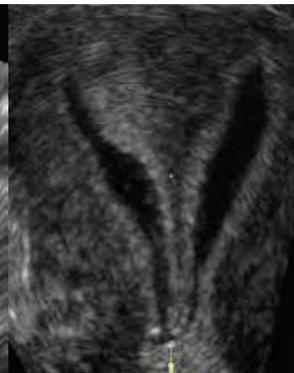

3D sonohysterogram of a complete uterine septum

Almost half of all endometrial polyps seen on SHG were missed using solely transvaginal US

Hartman A, Hartman, M. Sonohysterography versus Transvaginal Sonography. A blinded retrospective case series on 800 consecutive patients with polyps diagnosed on Sonohysterography. How accurate is Transvaginal Sonography? American Institute of Ultrasound in Medicine Conference, March 2010.

Sonohysterography is the method of choice to evaluate broids

Becker E Jr., Lev-Toa A.S, Kaufman EP, Halpern EJ, Edelweiss MI, Kurtz AB. The added value of transvaginal sonohysterography over transvaginal sonography alone in women with known or suspected leiomyoma. J Ultrasound Med. 2002 Mar; 21(3): 237-47.

Sonohysterography is an important part of the routine investigational process for evaluating infertility.

Tur-Kaspa I, Gal M, Hartman M, Hartman J, Hartman, A. A prospective evaluation of uterine abnormalities by saline infusion sonohysterography in 1009 women with infertility or abnormal uterine bleeding. Fertility and Sterility® December 2006; Vol. 86, No. 6.



ver 07/28/2016

BATRIK KRYSAL VIEW™ Ultrasound Gel

KRYSAL VIEW Ultrasound Gel is a conductive medium that creates a continuous bond between the skin and the probe transducer, eliminating air space and helping to produce a clear image. KRYSAL VIEW Ultrasound Gel is fluid enough to fill all empty spaces and viscous enough to stay in place, consistently delivering superior image quality.

- Safe on probes due to the high viscosity of the active agents and vacuum technology used in its production
- Offered in medium or high viscosity
- Non-greasy, water-soluble, non-staining formulation
- Unaffected by salt on the skin
- Hypoallergenic, paraben-free (CIT/MIT and paraben-free options available, which are preservative that some patients, might be allergic to)
- Non-toxic
- Odourless

BATRIK KRYSAL VIEW™ Probe Cleaning Solution

Removes gel residue on the probe for clearer imaging.

It also:

- Prevents cross-contamination
- Alcohol-free
- Extends the life of the probe



Item #	Description		Packaging
GL-ST0020	Ultrasound Gel Single-Use Pouch, Sterile	20 mL	240 Units / Box
GL-0020	Ultrasound Gel Single-Use Pouch, Non-Sterile	20 mL	240 Units / Box
GL-0250	Ultrasound Gel Bottle, Non-Sterile	250 mL	40 / Box
GL-1000	Ultrasound Gel Bottle, Non-Sterile	1,000 mL	12 / Box
GL-5000	Ultrasound Gel Cubitainer, Non-Sterile	5,000 mL	4 / Box
SS-22-0250	Ultrasound Probe Cleaning Solution Bottle	250 mL	25 / Box

ver 05/26/2020



BATRIK MEDICAL MANUFACTURING

850 Halpern Ave., Montreal (Dorval), QC H9P 1G6

P.: 514-631-7988 • T.-F.: 1-800-661-5432 • F.: 514-631-9083 • E.: info@surgmed.com • W.: www.surgmed.com